



**System provider
for Incontinence Care**

Introduction

With our assortment of revulsive incontinence products, we as GHC German Health Care GmbH focus on an area, where our mission is to offer our customers an extensive range of different products as well as offering a very high level of expertise. And our customers can rely on that – Even during challenging times!

Within our product assortment, you will find a wide selection of transurethral and suprapubic balloon catheters, different catheter insertion kits, a wide range of urine collection systems and bags as well as bag holding systems and accessories.

We have positioned ourselves well for the future by having expanded our warehouse capacities, added a new warehouse management software as well as keeping up the high level of service orientation with various service offers. For example, we have created a customer portal where all relevant MDR-documents are available for you to view or download at any time. This service saves you a lot of time and gives you the benefit of always having access to all necessary documents.

In our current catalogue you will find new products that we have added to our assortment for our customers ranging from homecare companies, medical supply stores and clinics. For instance, we are introducing our new CARE *CON* series (page 38-40), which are Male External Catheters (MEC) in a very good and reliable quality at an attractive price. Further, we are introducing our new transurethral soft silicone balloon catheter (page 9), which due to the softer shaft material is an excellent alternative to the current balloon catheter range in the homecare market. Alongside that, we are implementing an additional new transurethral silicone catheter with a flute tip (page 14). This variant is a good solution for patients with problematic urine findings. Lastly, we also have implemented new suprapubic balloon catheters as alternatives (pages 25 and 29) to our previous product line-up. Besides the high product quality, we also wanted to enhance our delivery capabilities by adding the new products to our assortment.

With this catalogue we offer you an overview of the complete range of GHC's product portfolio. You are warmly invited to visit **www.ghc-medical.de** for up-to-date information.

If our offer has sparked your interest, feel free to get in touch with us.

We look forward talking to you.

Warmest greetings from the GHC management,


Dirk Dohmann
Kay Kruse

Table of Contents

Introduction	page 3
Table of Contents	page 4,5
Quality policy of GHC	page 6
Catheter guide	page 7
CARE FLOW Transurethral catheter	
CARE FLOW 2-way 100% silicone catheter, nelaton tip	page 8
CARE FLOW SOFT- 2-way Soft 100% silicone catheter, nelaton tip new	page 9
CARE FLOW 2-way integral 100% silicone catheter, nelaton tip and plug new	page 10
BLUECATH 2-way integral 100% silicone catheter, nelaton tip and plug new	page 11
CARE FLOW 2-way silicone elastomer catheter, nelaton tip	page 12
CARE FLOW 2-way 100% silicone catheter, tiemann tip new	page 13
CARE FLOW 2-way 100% silicone catheter, flute tip new	page 14
CARE FLOW 2-way silicone elastomer catheter, couvelaire tip	page 15
CARE FLOW 2-way siliconised latex catheter, nelaton tip	page 16
CARE FLOW PVC single-use female catheter, nelaton tip	page 17
CARE FLOW PVC single-use male catheter, nelaton and tiemann tip	page 18
CARE FLOW Irrigation catheter	
CARE FLOW 3-way 100% silicone irrigation catheter, nelaton tip	page 19
CARE FLOW 3-way 100% silicone irrigation catheter, dufour tip	page 20
CARE FLOW 3-way 100% silicone irrigation catheter, couvelaire tip	page 21
CARE FLOW 3-way 100% silicone irrigation catheter, mercier tip	page 22
CARE FLOW Suprapubic catheter	
CARE FLOW UNIBAL 2-way integral 100% silicone catheter, central open tip	page 23
CARE FLOW UNIBAL 2-way integral 100% silicone catheter, central open tip, with plug	page 24
CARE FLOW 2-way 100% silicone catheter integral, central open tip, with plug new	page 25
CARE FLOW 2-way integral 100% silicone catheter, central open tip, short	page 26
BLUECATH 2-way integral 100% silicone catheter, dilatation tip, with plug	page 27
BLUECATH 2-way integral 100% silicone catheter, central open tip, with plug	page 28
CARE FLOW 2-way 100% silicone catheter integral blue, central open tip, with plug new	page 29
BLUECATH 2-way integral 100% silicone catheter, short, with plug	page 30
CARE FLOW Nephrostomy catheter	page 31
CARE FLOW 2-way integral 100% silicone catheter, with dilatation tip	page 32
CARE FLOW Nephrostomy catheter	
CARE FLOW 2-way integral 100% silicone catheter, with dilatation tip	page 33

Table of Contents

CARE SET Catheter insertion kit and additional equipment	
Catheter insertion kit Basic and Clinic	page 34
Catheter insertion kit Comfort and Comfort Plus	page 35
Catheter insertion kit Comfort Eco	page 36
Additional equipment	page 37
CARE /SC Intermittent Catheter	
CARE CATH <i>Leir</i> IC disposable catheter, hydrophilic, 17 cm new	Seite 38
CARE CATH <i>Leir</i> IC disposable catheter, hydrophilic, 41 cm new	Seite 49
CARE CON Male External Catheter	
CARE CON Male External Catheter comfort / type I new	page 40
CARE CON Male External Catheter special / type II new	page 41
CARE CON Male External Catheter active / type III new	page 42
CARE FLOW Urine drainage systems and accessories	
CARE FLOW UB 4000 TUR Closed TUR irrigation bag system	page 43
CARE FLOW UB 4000 Closed urine drainage system	page 44
CARE FLOW UB 2000 Closed urine drainage system, incl. CARE BAG UB 2000	page 45
CARE FLOW UB 2000A Sterile single-use bedside bag	page 46
CARE FLOW UB 2000E Sterile single-use bedside bag, incl. UB 2000US	page 47
CARE FLOW UB 2000USo, UB 1500USo, SB 1500USo	page 48
CARE FLOW UB 1000 Closed urine drainage system for day & night	page 49
CARE FLOW UB 750/750US 1-chamber leg bag system	page 50
CARE FLOW UB 750F/750UF 1-chamber leg bag systems with flip valve	page 51
CARE FLOW UB 500/500US 3-chamber leg bag system	page 52
CARE FIX Leg bag fixation	
CARE FIX Universal Velcro strap	page 53
CARE FIX Velcro strap premium set, with and without buttons	page 54
CARE FIX Leg pants with integrated leg bag pocket	page 55
CARE FIX Fabric leg bag holder for the thigh	page 56
CARE FIX Lower leg cuffs	page 57
CARE FIX Premium one leg pants	page 58
CARE FLOW Catheter valve	
CARE FLOW and CARE VENT	page 59
CARE DRESS Dressing material	
Suprapubic dressing set	page 60
Abdominal PEG dressing set	page 61

German Health Care GmbH stands for innovative solutions in medical care of persons with physical disabilities. As an importer and distributor of products especially for the bladder management, we offer clinics, homecare providers and pharmacies a broad portfolio and individual solutions to improve the quality of life of your customers – whilst giving you the added safety of receiving high product quality as well as reliability.

Our own, extensive product range is thanks to the strategic cooperation with well-known manufacturers supplemented by established products. Due to our continuous exchange with our customers and the end-users, our products are not only in a state of steady improvement, but also new products are added to our assortment. Hence not only market-leading homecare providers value our high-quality product and services since several years, but also clinics and the purchasing communities of medical trade suppliers.

The production of GHC-products is guided by clear specifications of our manufacturers, based at different locations worldwide. To ensure the quality of the production processes, every production site is getting audited regularly. The manufacturers conduct quality controls directly after the production and before the goods gets shipped. The final inspections are conducted when the goods arrive at our high-bay warehouse in Wees. Before the GHC products are shipped to our customers, they get checked and tested extensively. Only perfectly produced products are getting released and our customers can rely on this strict quality policy.



This strict quality policy is also embedded in our quality management system according to DIN EN ISO 13485.

Not only can you rely on the quality of our products, but rather also on our employees, as they are the basis of our success. Our qualified medical device consultants in the back office as well as our field staff will gladly advise you on all questions concerning the GHC product portfolio and take care of it to process your requests as quickly as possible. Through continuous in-house workshops and further education as well as IT-supported processes we are guaranteeing this expertise, so that your inquiries, your orders, and your goods are processed as quickly as possible. You can count on that.

Ensure your incontinence supply quality by cooperating with GHC German Health Care GmbH, your system provider for revulsive incontinence care.



Size & Color Coding: 1 Charrière (CH = 1 French 0.33mm)

We carry catheters of different materials with different specifications to meet all the needs of our customers. Among them are catheters made of full silicone and Silicone elastomer for long-term drainage, but also latex catheter for short-term drainage.

Classic catheter materials



Catheter tips for transurethral and suprapubic urine drainage




Pictures may differ from the original

CARE FLOW
Transurethral 100% silicone catheter
2-way balloon catheter with nelaton tip

Product specifications:

- ✓ In CH sizes 06 – 10 the catheter comes with a stabilising guiderail.
- ✓ Due to the high balloon capacity of up to 30 ml, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The X-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long-term drainage.

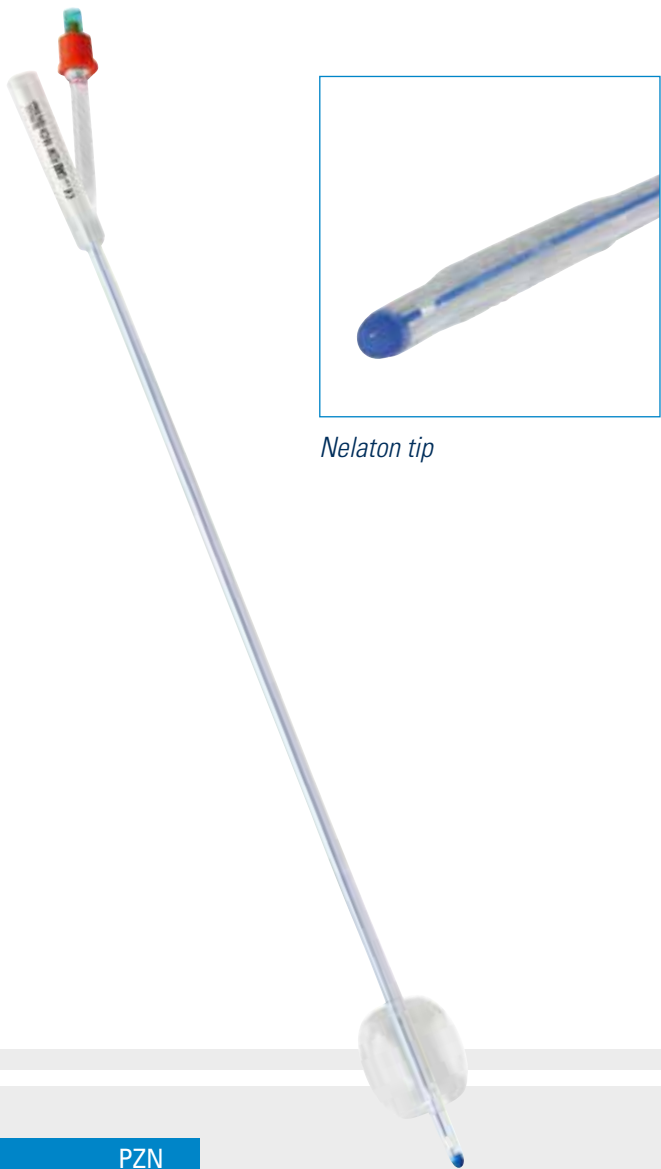

Nelaton tip



Balloon volume: 03-30 ml

Size	Balloon	Length	Article number	PZN
CH 06	03 ml	27 cm	01/01/11/06/03/GHC	01150196
CH 08	05 ml	27 cm	01/01/11/08/05/GHC	01150204
CH 10	05 ml	27 cm	01/01/11/10/05/GHC	01150210
CH 12	10 ml	40 cm	01/01/11/12/10/GHC	01150368
CH 14	10 ml	40 cm	01/01/11/14/10/GHC	01150374
CH 16	30 ml	40 cm	01/01/11/16/30/GHC	01150380
CH 18	30 ml	40 cm	01/01/11/18/30/GHC	01150397
CH 20	30 ml	40 cm	01/01/11/20/30/GHC	01150428
CH 22	30 ml	40 cm	01/01/11/22/30/GHC	01150434
CH 24	30 ml	40 cm	01/01/11/24/30/GHC	01150440
CH 26	30 ml	40 cm	01/01/11/26/30/GHC	01150457

Packaging Unit: 10 pieces
Packaging: sterile

Pictures may differ from the original

CARE FLOW
Transurethral *Soft* 100% silicone catheter
2-way balloon catheter with nelaton tip

Product specifications:

- ✓ In CH sizes 06 – 10 the catheter comes with a stabilising guiderail.
- ✓ Very high wearing comfort, due to the soft-silicone material of the catheter shaft.
- ✓ Safe fixation of the catheter in the bladder due to the balloon.
- ✓ Use our blocking syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The X-ray line allows for radiological imaging of the catheter.
- ✓ The catheter can be used for long-term drainage.



Nelaton tip



Balloon volume: 03-30 ml

Size	Balloon	Length	Article number	PZN
CH 06	03 ml	31 cm	01/01/01/06/03/GHC	17883185
CH 08	05 ml	31 cm	01/01/01/08/05/GHC	17883216
CH 10	05 ml	31 cm	01/01/01/10/05/GHC	17883222
CH 12	10 ml	40 cm	01/01/01/12/10/GHC	17883239
CH 14	10 ml	40 cm	01/01/01/14/10/GHC	17883251
CH 16	30 ml	40 cm	01/01/01/16/30/GHC	17883268
CH 18	30 ml	40 cm	01/01/01/18/30/GHC	17883274
CH 20	30 ml	40 cm	01/01/01/20/30/GHC	17883280
CH 22	30 ml	40 cm	01/01/01/22/30/GHC	17883297
CH 24	30 ml	40 cm	01/01/01/24/30/GHC	17883305
CH 26	30 ml	40 cm	01/01/01/26/30/GHC	17883311

Packaging Unit: 10 pieces
Packaging: sterile

Pictures may differ from the original

CARE FLOW

Transurethral *integral* 100% silicone catheter

2-way integral balloon catheter with nelaton tip and plug

- Product specifications:
- ✓ Safe fixation of the catheter, due to the seamlessly incorporated *integral* balloon into the shaft.
 - ✓ Use our blocker syringe containing 10% glycerine solution here.
 - ✓ Excellent biocompatibility and high resistance to encrustation.
 - ✓ The X-ray line allows for radiological imaging of the catheter.
 - ✓ MRT – compatibility.
 - ✓ The catheter can be used for long-term drainage.



Nelaton tip,
integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	41 cm	01/01/25/12/05/GHC	18064437
CH 14	05 ml	41 cm	01/01/25/14/05/GHC	18064443
CH 16	10 ml	41 cm	01/01/25/16/10/GHC	18064466
CH 18	10 ml	41 cm	01/01/25/18/10/GHC	18064472
CH 20	10 ml	41 cm	01/01/25/20/10/GHC	18064489
CH 22	10 ml	41 cm	01/01/25/22/10/GHC	18064495

Packaging Unit: 10 pieces

Packaging: sterile




Balloon volume: 05-10 ml

BLUECATH

Transurethral *integral* 100% silicone catheter

2-way integral balloon catheter with nelaton tip and plug

- Product specifications:
- ✓ Safe fixation of the catheter, due to the seamlessly incorporated *integral* balloon into the shaft.
 - ✓ Use our blocker syringe containing 10% glycerine solution here.
 - ✓ Excellent biocompatibility and high resistance to encrustation.
 - ✓ MRT – compatibility.
 - ✓ The catheter can be used for long-term drainage.

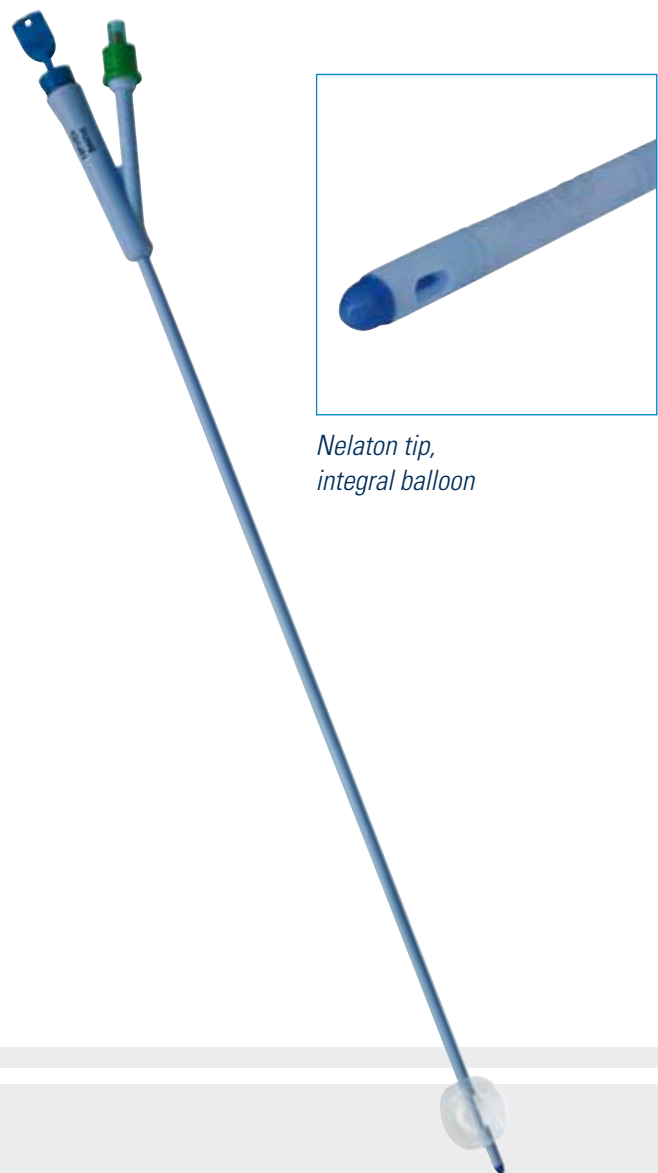

Nelaton tip,
integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	41 cm	01/01/26/12/05/GHC	18064503
CH 14	05 ml	41 cm	01/01/26/14/05/GHC	18064526
CH 16	10 ml	41 cm	01/01/26/16/10/GHC	18064532
CH 18	10 ml	41 cm	01/01/26/18/10/GHC	18064549
CH 20	10 ml	41 cm	01/01/26/20/10/GHC	18064555
CH 22	10 ml	41 cm	01/01/26/22/10/GHC	18064561

Packaging Unit: 10 pieces

Packaging: sterile



Balloon volume: 05-10 ml


CARE FLOW

Transurethral silicone elastomer catheter

2-way balloon catheter with nelaton tip

Product specifications:

- ✓ High elasticity and wearing comfort, due to the silicone coating on the inner and outer catheter wall.
- ✓ Due to the high balloon capacity, the catheter is safety fixed in the bladder.
- ✓ For blocking use our blocker syringe containing 10% glycerine solution.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The catheter can be used for long-term drainage.

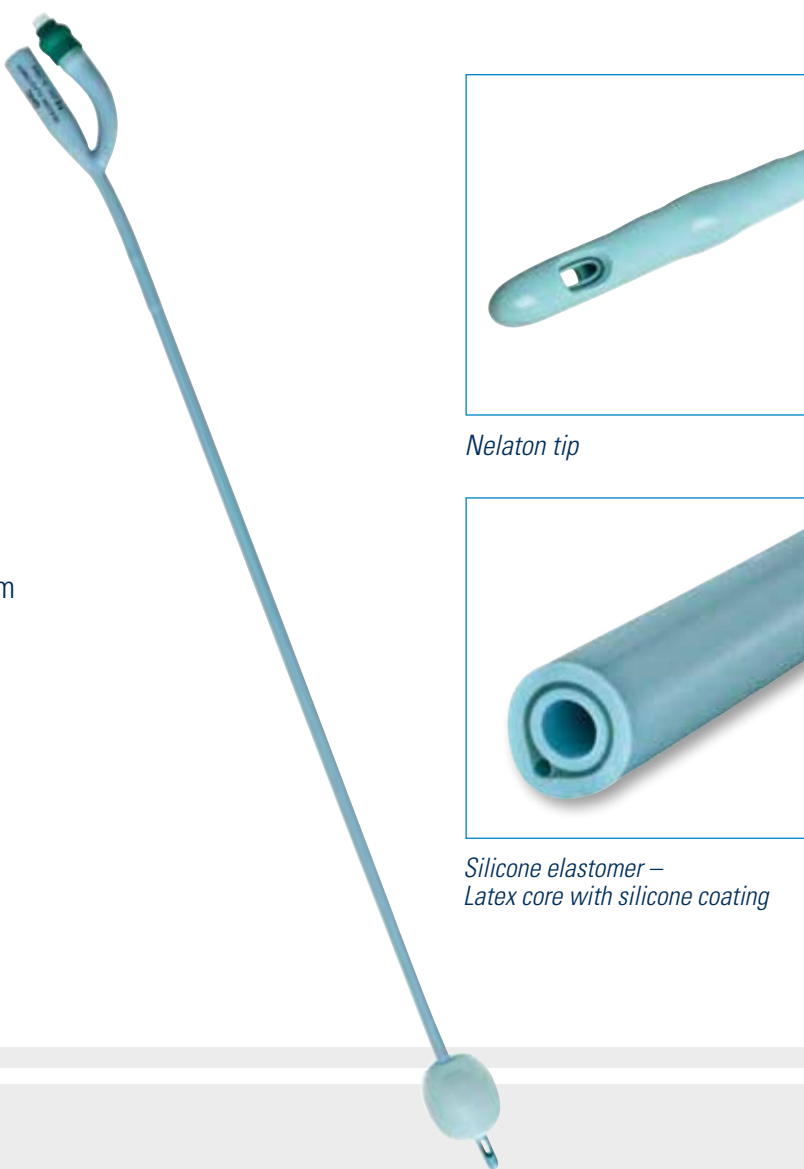

Nelaton tip



Silicone elastomer – Latex core with silicone coating



Balloon volume: 10-30 ml

Size	Balloon	Length	Article number	PZN
CH 12	10 ml	40 cm	01/01/10/12/10/GHC	01149371
CH 14	30 ml	40 cm	01/01/10/14/30/GHC	01149388
CH 16	30 ml	40 cm	01/01/10/16/30/GHC	01149394
CH 18	30 ml	40 cm	01/01/10/18/30/GHC	01149939
CH 20	30 ml	40 cm	01/01/10/20/30/GHC	01149945
CH 22	30 ml	40 cm	01/01/10/22/30/GHC	01149951
CH 24	30 ml	40 cm	01/01/10/24/30/GHC	01149968
CH 26	30 ml	40 cm	01/01/10/26/30/GHC	01150173

Packaging Unit: 10 pieces
Packaging: sterile

Pictures may differ from the original


CARE FLOW

Transurethral 100% silicone catheter

2-way balloon catheter with tiemann tip

Product specifications:

- ✓ An atraumatic dilating ball tip ensures a smooth passage of the prostatic fossa.
- ✓ Safe fixation of the catheter due to its high balloon capacity of up to 30 ml.
- ✓ Use our blocking syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The X-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility
- ✓ The catheter can be used for long term drainage.

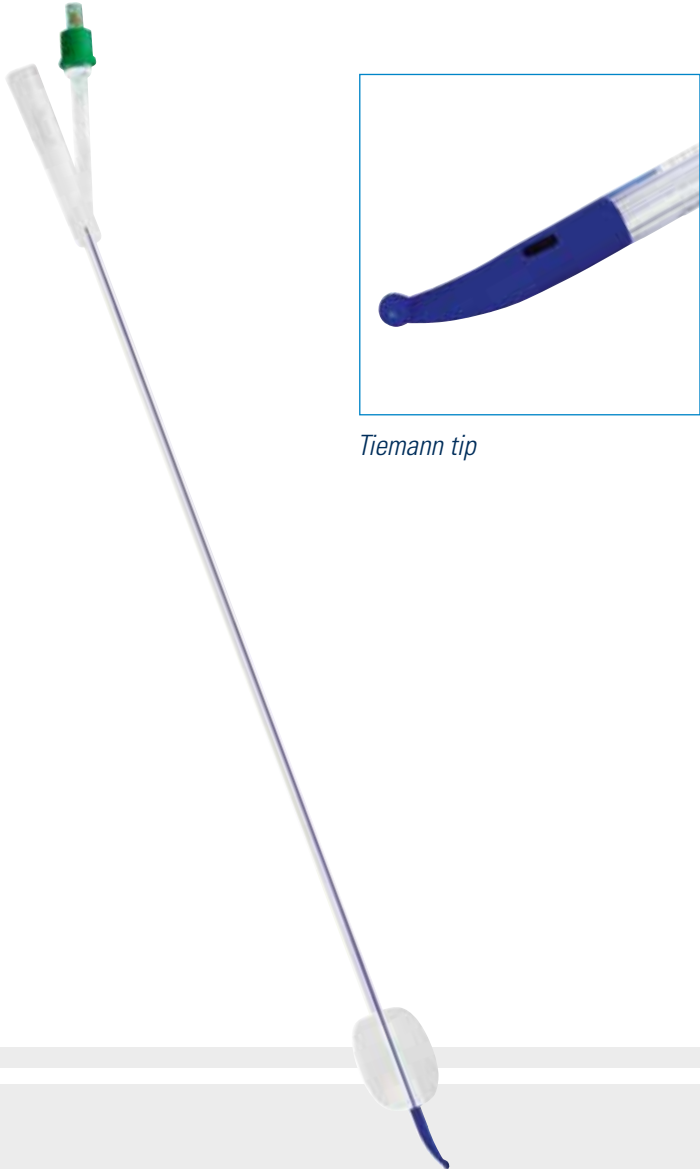

Tiemann tip



Balloon volume: 5-30 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/01/08/12/05/GHC	17933670
CH 14	20 ml	40 cm	01/01/08/14/20/GHC	17933687
CH 16	30 ml	40 cm	01/01/08/16/30/GHC	17933693
CH 18	30 ml	40 cm	01/01/08/18/30/GHC	17933701
CH 20	30 ml	40 cm	01/01/08/20/30/GHC	17933718
CH 22	30 ml	40 cm	01/01/08/22/30/GHC	17933724
CH 24	30 ml	40 cm	01/01/08/24/30/GHC	17933730

Packaging Unit: 10 pieces
Packaging: sterile

Pictures may differ from the original


CARE FLOW

Transurethral 100% silicone catheter

2-way balloon catheter with flute tip

Product specifications:

- ✓ Excellent drainage performance, due to the central open tip, drainage eye, and high lumen, even for problematic urine findings.
- ✓ An atraumatic flute tip ensures an easy placement of the catheter.
- ✓ Safe fixation of the catheter due to its high balloon capacity of up to 30 ml.
- ✓ Use our blocking syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The X-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long term drainage.

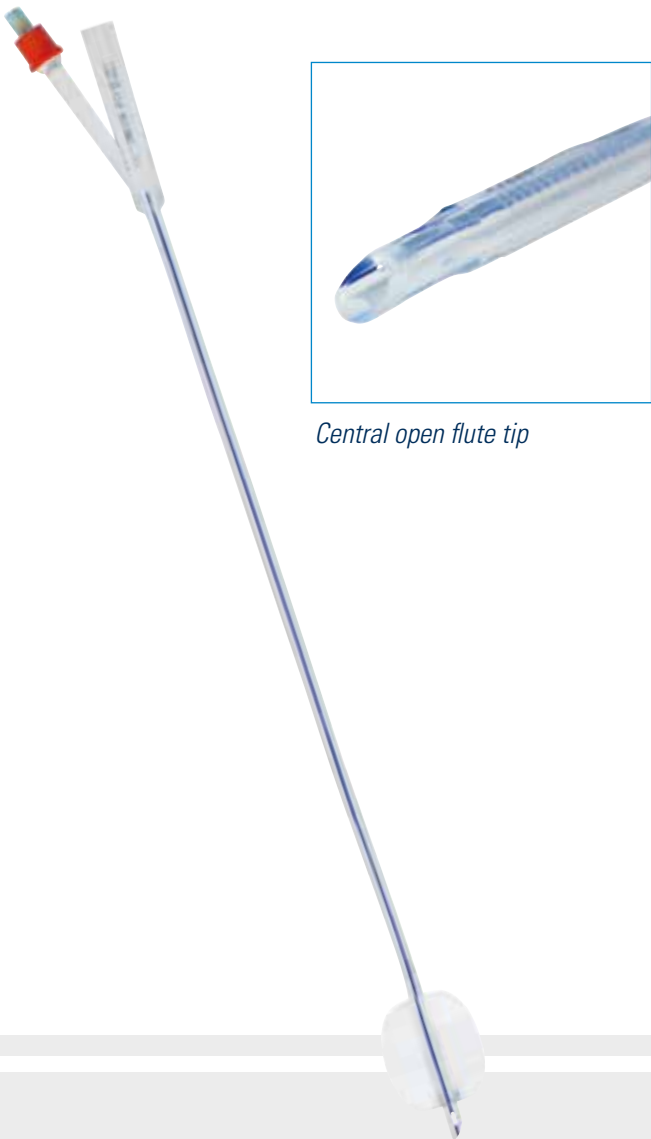

Central open flute tip

Size	Balloon	Length	Article number	PZN
CH 12	10 ml	40 cm	01/01/07/12/10/GHC	17933747
CH 14	10 ml	40 cm	01/01/07/14/10/GHC	17933753
CH 16	30 ml	40 cm	01/01/07/16/30/GHC	17933776
CH 18	30 ml	40 cm	01/01/07/18/30/GHC	17933782
CH 20	30 ml	40 cm	01/01/07/20/30/GHC	17933799
CH 22	30 ml	40 cm	01/01/07/22/30/GHC	17933807

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 10-30 ml


CARE FLOW

Transurethral silicone elastomer catheter


2-way balloon catheter with couvelaire tip

Product specifications:

- ✓ High elasticity and wearing comfort, due to the silicone coating on the inner and outer catheter wall.
- ✓ Optimal drainage and irrigation due to two staggered drainage eyes and open couvelaire tip.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Safe fixation of the catheter, due to the high balloon capacity of up to 60 ml.
- ✓ For blocking use our blocker syringe containing 10% glycerine solution.
- ✓ The catheter can be used for long-term drainage.

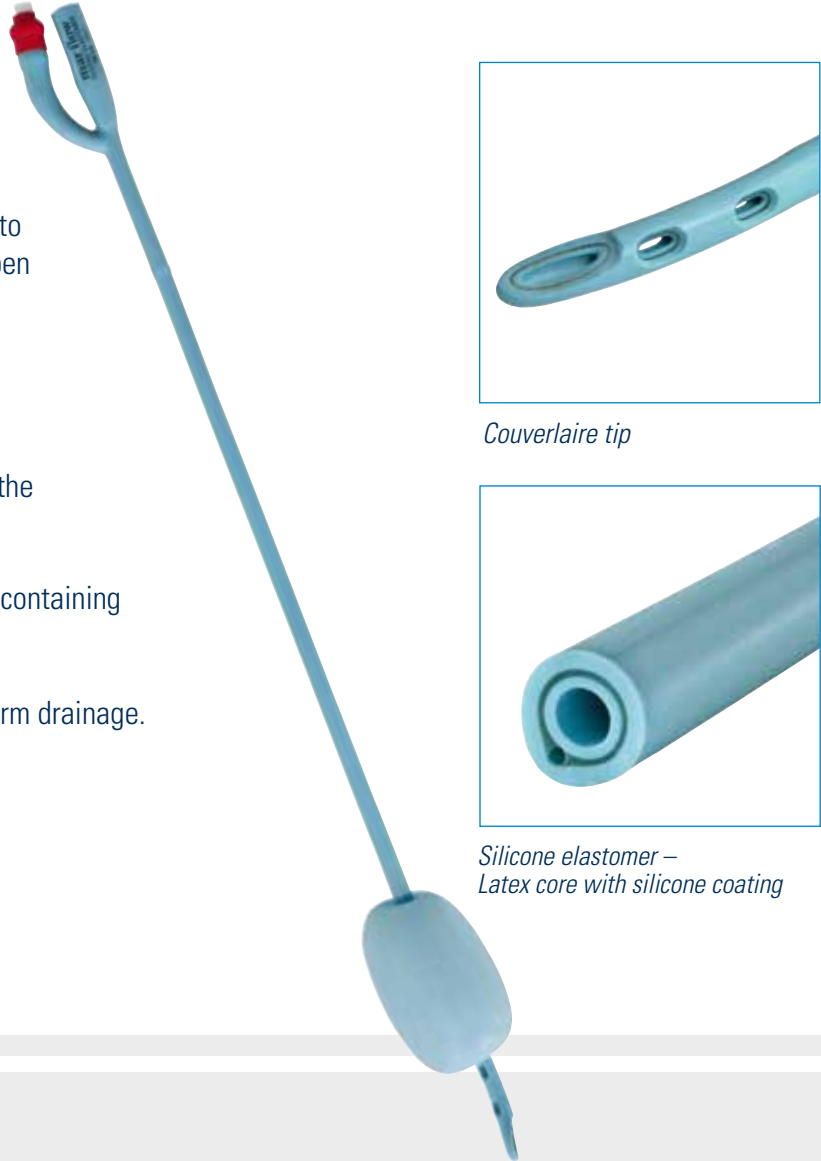

Couvelaire tip

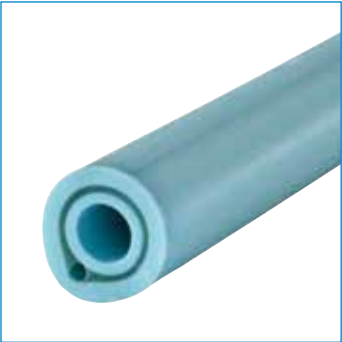

Silicone elastomer – Latex core with silicone coating

Size	Balloon	Length	Article number	PZN
CH 18	60 ml	40 cm	01/01/15/18/60/GHC	06434975
CH 20	60 ml	40 cm	01/01/15/20/60/GHC	06434952
CH 22	60 ml	40 cm	01/01/15/22/60/GHC	06434969
CH 24	60 ml	40 cm	01/01/15/24/60/GHC	09710080

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: up to 60 ml


CARE FLOW

Transurethral siliconised latex catheter

2-way balloon catheter with nelaton tip

Product specifications:

- ✓ Made of high-quality, siliconized medical latex.
- ✓ The catheter offers high patient comfort.
- ✓ Due to the high balloon capacity, the catheter is safety fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The catheter is suited for short-term drainage of up to 5 days only.
- ✓ This product contains natural latex which may cause allergic reactions.

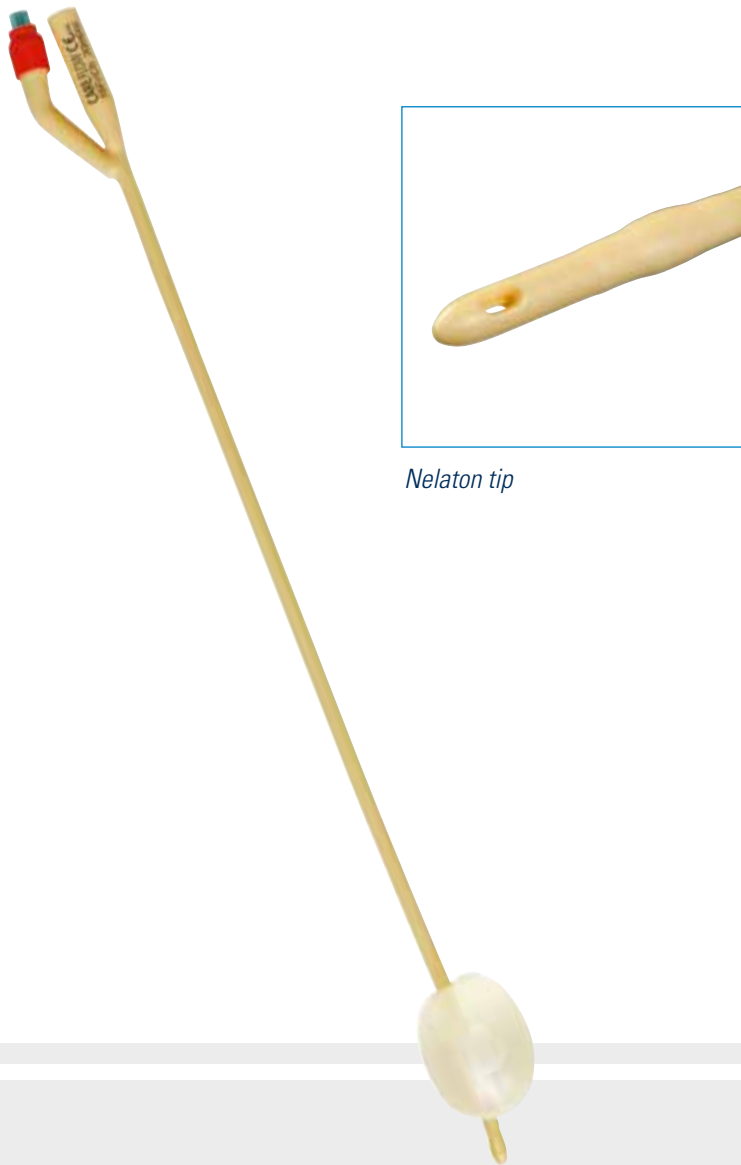

Nelaton tip

Size	Balloon	Length	Article number	PZN
CH 12	10 ml	40 cm	01/01/13/12/10/GHC	01150701
CH 14	10 ml	40 cm	01/01/13/14/10/GHC	01150859
CH 16	30 ml	40 cm	01/01/13/16/30/GHC	01150865
CH 18	30 ml	40 cm	01/01/13/18/30/GHC	01150871
CH 20	30 ml	40 cm	01/01/13/20/30/GHC	01150888
CH 22	30 ml	40 cm	01/01/13/22/30/GHC	01150919

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 10-30 ml

Pictures may differ from the original


CARE FLOW

Transurethral PVC single-use female catheter

Transparent with funnel and nelaton tip, DEHP-free

Product specifications:

- ✓ Disposable female catheter with two staggered atraumatic drainage eyes.
- ✓ Rounded and closed tip.
- ✓ Soft PVC material.
- ✓ DEHP-free.
- ✓ The catheter is for single use only.

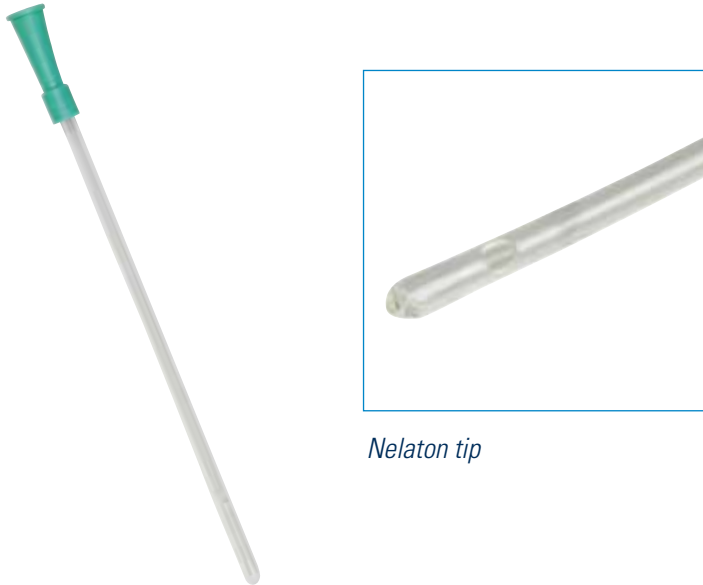

Nelaton tip

Size	Length	Article number	PZN
CH 08	20 cm	01/01/35/08/00/GHC	18320975
CH 10	20 cm	01/01/35/10/00/GHC	18320981
CH 12	20 cm	01/01/35/12/00/GHC	18320998
CH 14	20 cm	01/01/35/14/00/GHC	18321006
CH 16	20 cm	01/01/35/16/00/GHC	18321035

Packaging Unit: 100 pieces
Packaging: sterile



Laterally displaced eyes

Pictures may differ from the original

CARE FLOW

Transurethral PVC single-use male catheter

Transparent with funnel, nelaton and tiemann tip

Product specifications:

- ✓ Disposable male catheter with two staggered atraumatic drainage eyes.
- ✓ Rounded and closed Nelaton or Tiemann tip, for a smooth passage of the prostatic fossa.
- ✓ Soft PVC material.
- ✓ The catheter is for single use only.



Nelaton tip

Size	Length	Article number	PZN
CH 10	40 cm	01/01/20/10/00/GHC	08871349
CH 12	40 cm	01/01/20/12/00/GHC	08871355
CH 14	40 cm	01/01/20/14/00/GHC	08871361
CH 16	40 cm	01/01/20/16/00/GHC	08871378
CH 18	40 cm	01/01/20/18/00/GHC	08871384

Tiemann tip

Size	Length	Article number	PZN
CH 12	40 cm	01/01/21/12/00/GHC	01929620
CH 14	40 cm	01/01/21/14/00/GHC	01929637
CH 16	40 cm	01/01/21/16/00/GHC	01929666
CH 18	40 cm	01/01/21/18/00/GHC	01971792
CH 20	40 cm	01/01/21/20/00/GHC	09219740

Packaging Unit: CH 08-18: 50 pieces each, CH 20: 30 pieces
Packaging: sterile

Pictures may differ from the original

CARE FLOW

Transurethral 100% silicone irrigation catheter

3-way balloon catheter with nelaton tip

Product specifications:

- ✓ An atraumatic Nelaton tip ensures a smooth passage of the prostatic fossa.
- ✓ Optimal drainage and irrigation performance, due to two staggered drainage eyes.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for post-operative long-term drainage.



Size	Balloon	Length	Article number	PZN
CH 18	30 ml	40 cm	01/01/22/18/30/GHC	15618237
CH 20	30 ml	40 cm	01/01/22/20/30/GHC	15618243
CH 22	30 ml	40 cm	01/01/22/22/30/GHC	15618266
CH 24	30 ml	40 cm	01/01/22/24/30/GHC	15618295

Packaging Unit: 10 pieces
Packaging: sterile

Pictures may differ from the original



Balloon volume: 30 ml

CARE FLOW
Transurethral 100% silicone irrigation catheter
3-way balloon catheter with dufour tip

Product specifications:

- ✓ An atraumatic and kink resistant Defour tip ensures a smooth passage of the prostatic fossa.
- ✓ Optimal drainage and irrigation performance, due to two staggered drainage eyes and an open tip.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray tip allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for post-operative long-term drainage.



Dufour tip

Size	Balloon	Length	Article number	PZN
CH 18	40 ml	40 cm	01/01/16/18/40/GHC	01151178
CH 20	50 ml	40 cm	01/01/16/20/50/GHC	01151184
CH 22	60 ml	40 cm	01/01/16/22/60/GHC	01151190
CH 24	70 ml	40 cm	01/01/16/24/70/GHC	01151379

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 40-70 ml

Pictures may differ from the original

CARE FLOW
Transurethral 100% silicone irrigation catheter
3-way balloon catheter with couvelaire tip

Product specifications:

- ✓ Optimal drainage and irrigation due to two staggered drainage eyes and an open tip.
- ✓ An atraumatic and kink resistant Couvelaire tip ensures a smooth passage of the prostatic fossa.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray tip allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for post-operative long-term drainage.



Couvelaire tip

Size	Balloon	Length	Article number	PZN
CH 18	40 ml	40 cm	01/01/17/18/40/GHC	01151385
CH 20	50 ml	40 cm	01/01/17/20/50/GHC	01151391
CH 22	60 ml	40 cm	01/01/17/22/60/GHC	01151416
CH 24	70 ml	40 cm	01/01/17/24/70/GHC	01151439

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 40-70 ml

Pictures may differ from the original

CARE FLOW

Transurethral 100% silicone irrigation catheter

3-way balloon catheter with mercier tip

- Product specifications:
- ✓ An atraumatic and kink resistant Mercier tip ensures a smooth passage of the prostatic fossa.
 - ✓ Optimal drainage and irrigation due to two staggered drainage eyes and an open tip.
 - ✓ Excellent biocompatibility and high resistance to encrustation.
 - ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
 - ✓ Use our blocker syringe containing 10% glycerine solution here.
 - ✓ The x-ray tip allows for radiological imaging of the catheter.
 - ✓ MRT – compatibility.
 - ✓ The catheter can be used for post-operative long-term drainage.



Mercier tip

Size	Balloon	Length	Article number	PZN
CH 18	40 ml	40 cm	01/01/18/18/40/GHC	01151511
CH 20	50 ml	40 cm	01/01/18/20/50/GHC	01151528
CH 22	60 ml	40 cm	01/01/18/22/60/GHC	01151534
CH 24	70 ml	40 cm	01/01/18/24/70/GHC	01151540

Packaging Unit: 10 pieces

Packaging: sterile

Pictures may differ from the original

CARE FLOW UNIBAL

Suprapubic *integral* 100% silicone catheter

2-way integral balloon catheter

- Product specifications:
- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
 - ✓ The catheter offers high patient comfort.
 - ✓ High drainage performance, due to the high lumen and the central open tip.
 - ✓ We recommend using or blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter
 - ✓ Excellent biocompatibility and high resistance to encrustation.
 - ✓ Short and central open tip.
 - ✓ The x-ray line allows for radiological imaging of the catheter.
 - ✓ MRT – compatibility.



Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/12/10/03/GHC	01152568
CH 12	05 ml	40 cm	01/02/12/12/05/GHC	01152574
CH 14	05 ml	40 cm	01/02/12/14/05/GHC	01152605
CH 16	10 ml	40 cm	01/02/12/16/10/GHC	01152611
CH 18	10 ml	40 cm	01/02/12/18/10/GHC	01152628
CH 20	10 ml	40 cm	01/02/12/20/10/GHC	01152700
CH 22	10 ml	40 cm	01/02/12/22/10/GHC	01152717

Packaging Unit: 10 pieces

Packaging: sterile

Pictures may differ from the original

CARE FLOW UNIBAL

Suprapubic *integral* 100% silicone catheter

2-way integral balloon catheter with plug

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers high patient comfort.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ We recommend using or blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.



Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/14/10/03/GHC	10198115
CH 12	05 ml	40 cm	01/02/14/12/05/GHC	10198121
CH 14	05 ml	40 cm	01/02/14/14/05/GHC	10198138
CH 16	10 ml	40 cm	01/02/14/16/10/GHC	10198144
CH 18	10 ml	40 cm	01/02/14/18/10/GHC	10198150
CH 20	10 ml	40 cm	01/02/14/20/10/GHC	10198167
CH 22	10 ml	40 cm	01/02/14/22/10/GHC	10198173

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 03-10 ml

CARE FLOW

Suprapubic 100% silicone catheter *integral*

2-way balloon catheter with plug

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers an excellent wearing comfort.
- ✓ We recommend using our blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter
- ✓ 100% silicone for excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.



Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	41 cm	01/02/09/12/05/GHC	16840937
CH 14	05 ml	41 cm	01/02/09/14/05/GHC	16840943
CH 16	10 ml	41 cm	01/02/09/16/10/GHC	16840972
CH 18	10 ml	41 cm	01/02/09/18/10/GHC	16840989
CH 20	10 ml	41 cm	01/02/09/20/10/GHC	16840995
CH 22	10 ml	41 cm	01/02/09/22/10/GHC	16841003

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 05-10 ml

CARE FLOW UNIBAL

Suprapubic *integral* 100% silicone catheter

2-way integral balloon catheter, short

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Safe fixation of the catheter, due to the integral balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.



Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	33 cm	01/02/13/10/03/GHC	01152746
CH 12	05 ml	33 cm	01/02/13/12/05/GHC	01152812
CH 14	05 ml	33 cm	01/02/13/14/05/GHC	01152829
CH 16	10 ml	33 cm	01/02/13/16/10/GHC	01152841
CH 18	10 ml	33 cm	01/02/13/18/10/GHC	01152918
CH 20	10 ml	33 cm	01/02/13/20/10/GHC	01152924
CH 22	10 ml	33 cm	01/02/13/22/10/GHC	01152947

Packaging Unit: 10 pieces

Packaging: sterile



Balloon volume: 03-10 ml


BLUECATH

Suprapubic *integral* 100% silicone catheter

2-way integral balloon catheter with dilatation tip and plug

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ Safe fixation of the catheter, due to the integral balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Dilatation tip with central bore.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.

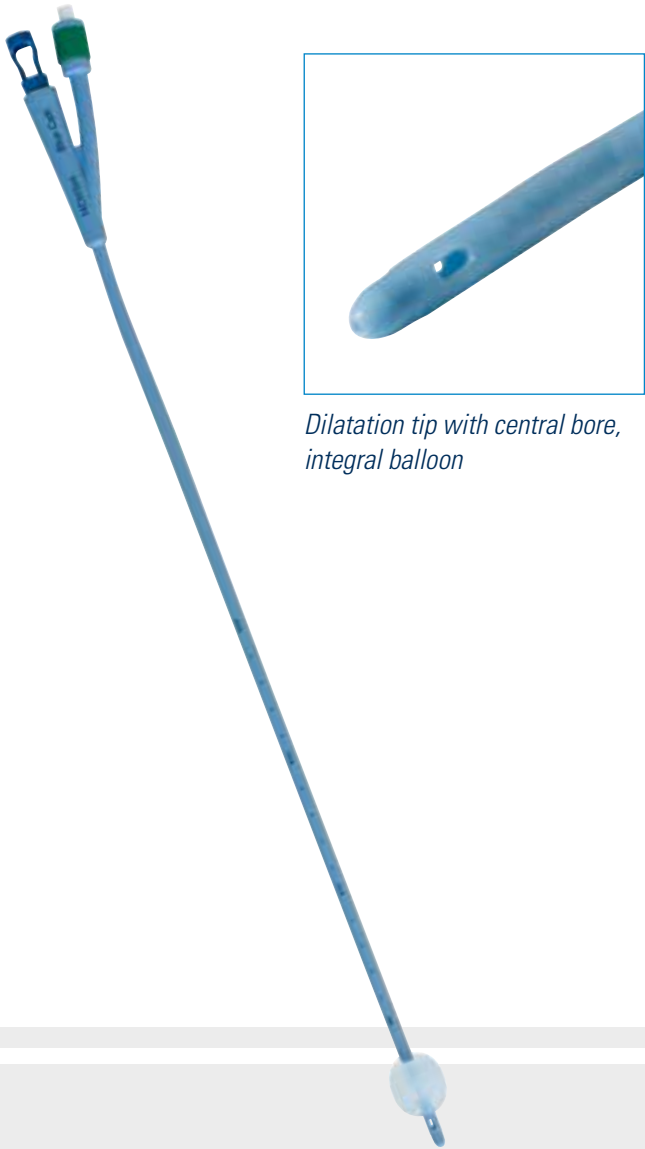

Dilatation tip with central bore, integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/02/21/12/05/GHC	11850948
CH 14	05 ml	40 cm	01/02/21/14/05/GHC	11850954
CH 16	10 ml	40 cm	01/02/21/16/10/GHC	11850983
CH 18	10 ml	40 cm	01/02/21/18/10/GHC	11851008
CH 20	10 ml	40 cm	01/02/21/20/10/GHC	11851014
CH 22	10 ml	40 cm	01/02/21/22/10/GHC	11851020

Packaging Unit: 10 pieces

Packaging: sterile




Balloon volume: 05-10 ml

BLUECATH
Suprapubic *integral* 100% silicone catheter
2-way integral balloon catheter with plug

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ Safe fixation of the catheter, due to the integral balloon.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short, central open tip.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.

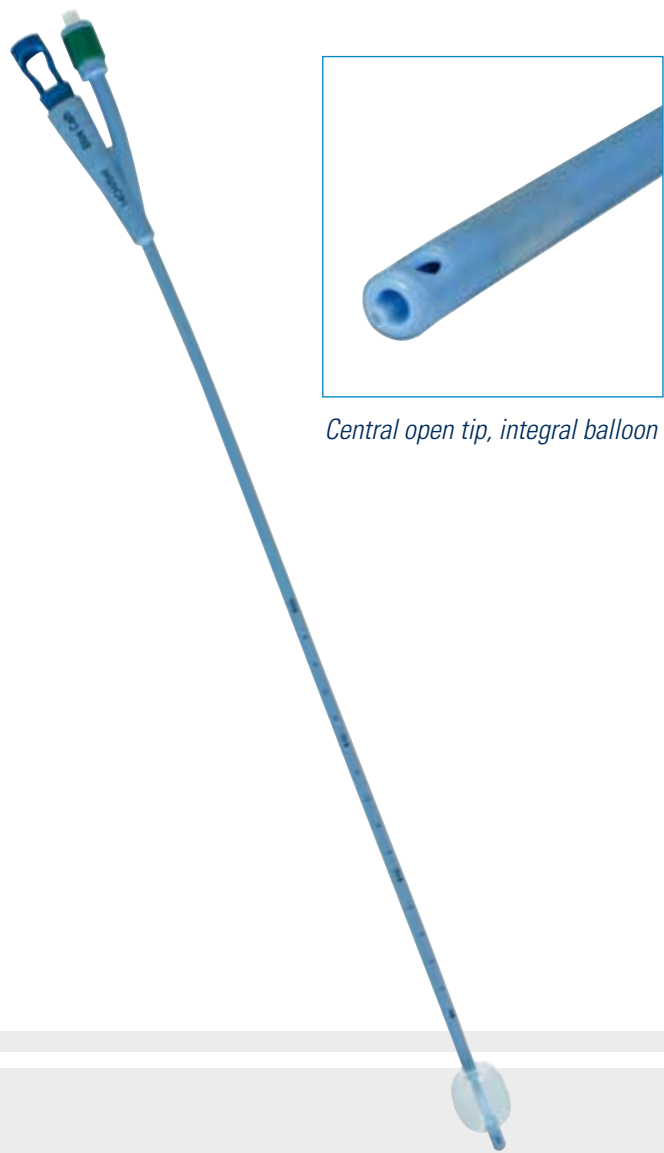

Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/02/22/12/05/GHC	11851037
CH 14	05 ml	40 cm	01/02/22/14/05/GHC	11851043
CH 16	10 ml	40 cm	01/02/22/16/10/GHC	11851072
CH 18	10 ml	40 cm	01/02/22/18/10/GHC	11851089
CH 20	10 ml	40 cm	01/02/22/20/10/GHC	11851095
CH 22	10 ml	40 cm	01/02/22/22/10/GHC	11851126
CH 24	10 ml	40 cm	01/02/22/24/10/GHC	11851132

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 05-10 ml

Pictures may differ from the original

CARE FLOW
Suprapubic 100% silicone catheter *integral blue*
2-way balloon catheter with plug

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers an excellent wearing comfort.
- ✓ We recommend using our blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter.
- ✓ 100% silicone for excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optimal placement.



Central open, integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	41 cm	01/02/11/12/05/GHC	16841026
CH 14	05 ml	41 cm	01/02/11/14/05/GHC	16841032
CH 16	10 ml	41 cm	01/02/11/16/10/GHC	16841049
CH 18	10 ml	41 cm	01/02/11/18/10/GHC	16841061
CH 20	10 ml	41 cm	01/02/11/20/10/GHC	16841078
CH 22	10 ml	41 cm	01/02/11/22/10/GHC	16841084

Packaging Unit: 10 pieces
Packaging: sterile



Balloon volume: 05-10 ml

Pictures may differ from the original

BLUECATH **Suprapubic *integral* 100% silicone catheter** 2-way integral balloon catheter with plug, short

Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short, central open tip.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.



Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	33 cm	01/02/23/12/05/GHC	11851149
CH 14	05 ml	33 cm	01/02/23/14/05/GHC	11851155
CH 16	10 ml	33 cm	01/02/23/16/10/GHC	11851161
CH 18	10 ml	33 cm	01/02/23/18/10/GHC	11851178
CH 20	10 ml	33 cm	01/02/23/20/10/GHC	11851184
CH 22	10 ml	33 cm	01/02/23/22/10/GHC	11851190

Packaging Unit: 10 pieces
Packaging: sterile




Balloon volume: 05-10 ml

CARE FLOW UNIBAL **Suprapubic puncture set** 2-way *integral* balloon catheter

Product specifications:

The puncture set contains:

- 100% silicone integral balloon catheter
- Stainless steel cannula with special cut
- Scalpel with small handle
- ✓ *Integral* balloon catheter with seamlessly incorporated balloon into the shaft.
- ✓ Optimal use of the inside diameter of the cannula due to integral balloon catheter: Charrière 10 and 12 with split cannula CH 13.8 (4.6 mm), puncture length 120 mm. Charrière 14 with split cannula CH 16.8 (5.6 mm), puncture length 120 mm.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Safe fixation of the catheter, due to the integral balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.

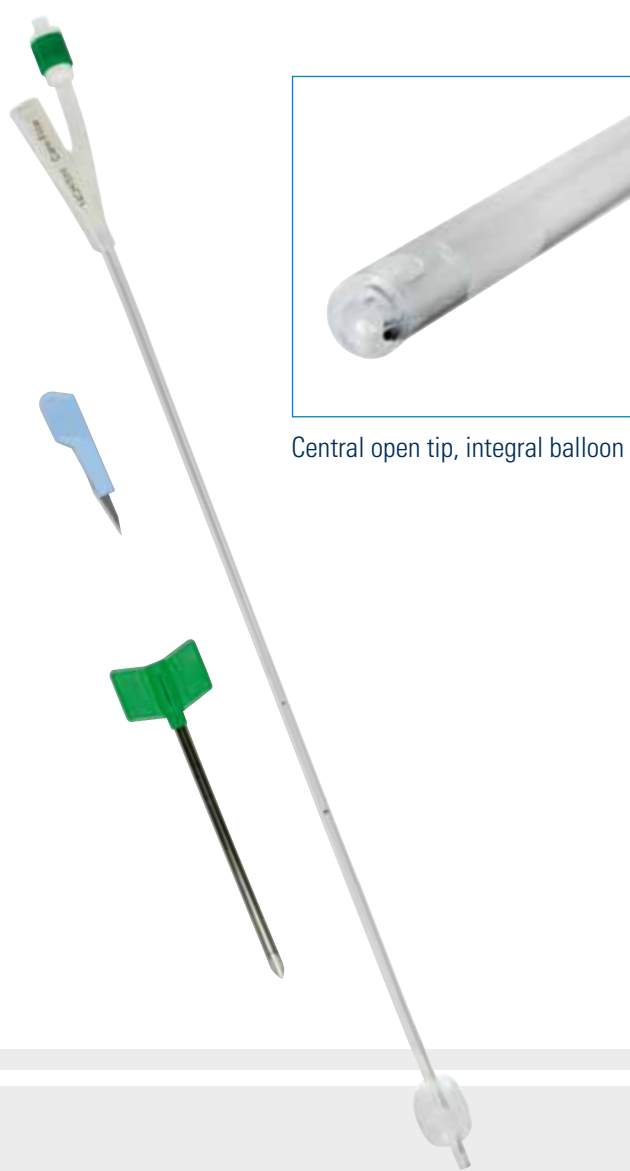

Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/15/10/03/GHC	01153059
CH 12	05 ml	40 cm	01/02/15/12/05/GHC	01153065
CH 14	05 ml	40 cm	01/02/15/14/05/GHC	01153088

Packaging Unit: 5 pieces
Packaging: sterile



Balloon volume: 03-05 ml


CARE FLOW UNIBAL

Suprapubic change set, incl. guidewire

2-way integral balloon catheter

Product specifications:

- This charge set contains:
- 100% silicone *integral* balloon catheter
 - Teflon-coated guidewire
- ✓ Integral balloon catheter with seamlessly incorporated balloon into the shaft.
 - ✓ High drainage performance, due to the high lumen and the central open tip.
 - ✓ Stable teflon-coated guidewire.
 - ✓ Safe fixation of the catheter, due to the integral balloon.
 - ✓ Use our blocker syringe containing 10% glycerine solution here.
 - ✓ The x-ray line allows for radiological imaging of the catheter.
 - ✓ MRT – compatibility.
 - ✓ Teflon-coated guidewire, lenght 80 cm, 3 cm straight and flexible tip, thickness 0,035".

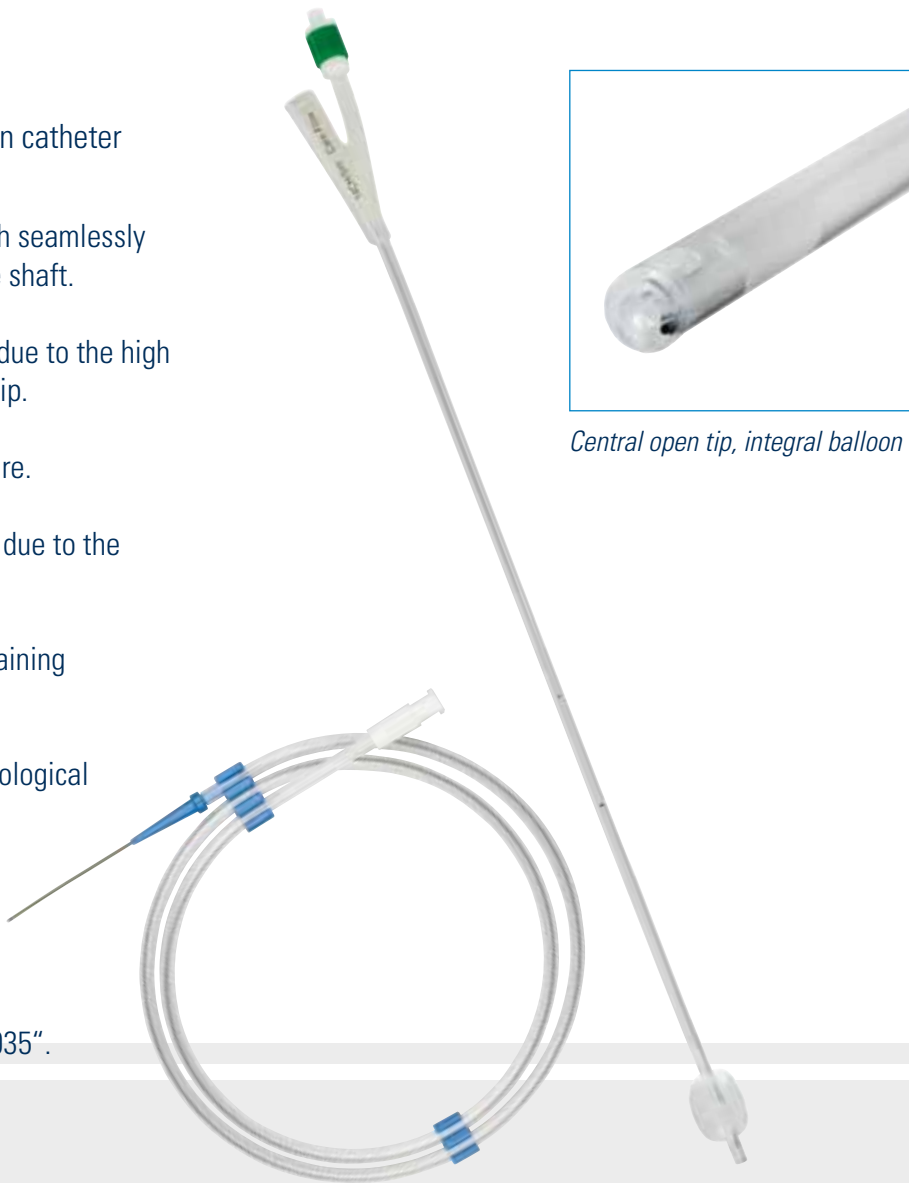

Central open tip, integral balloon

Guidewire

Wire	Tip	Length	Article number	PZN
coated	3 cm flex.	80 cm	01/02/18/80/35/GHC	14166431

Packaging Unit: 1 piece

Packaging: sterile

Change set, incl. Guidewire

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/16/10/03/GHC	01153154
CH 12	05 ml	40 cm	01/02/16/12/05/GHC	01153160
CH 14	05 ml	40 cm	01/02/16/14/05/GHC	01153183

Packaging Unit: 10 pieces

Packaging: sterile



Balloon volume: 03-05 ml


NEPHROSTOMY

Nephrostomy *integral* 100% silicone catheter

2-way integral balloon catheter with dilatation tip

Product specifications:

- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The catheters are suitable for long-term drainage.
- ✓ Dilating nelaton tip with central bore.
- ✓ High drainage performance, due to four drainage eyes located in front of and behind the balloon.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optimal placement.
- ✓ 100% silicone.

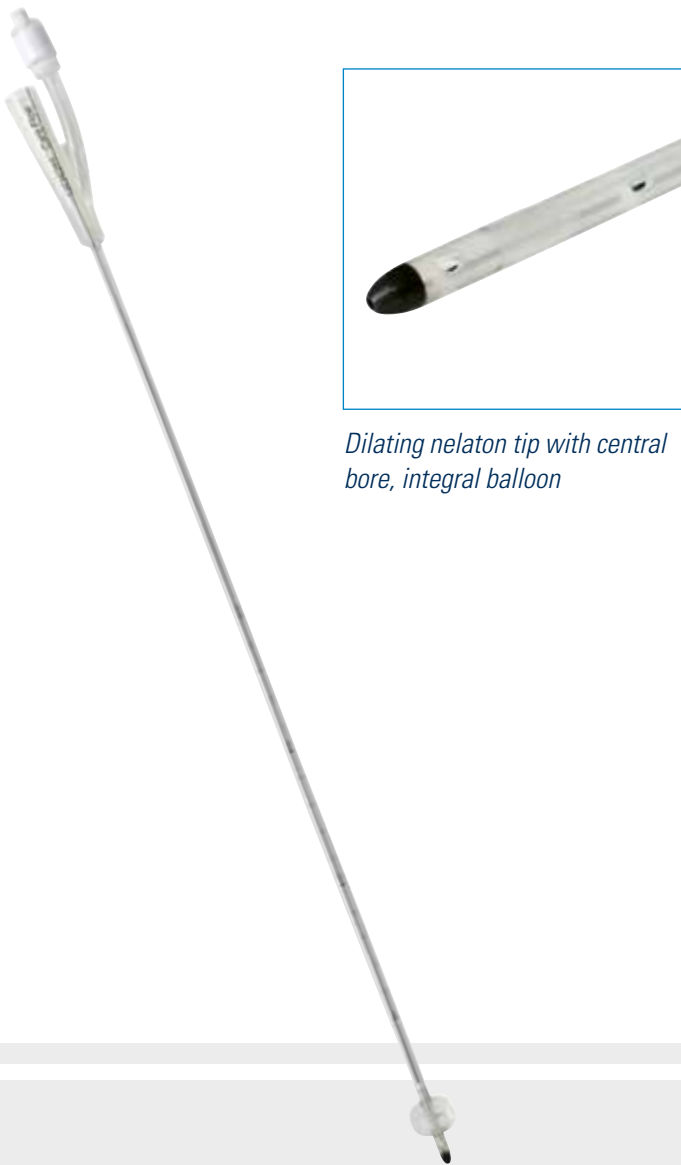


Dilating nelaton tip with central bore, integral balloon

Size	Balloon	Length	Article number	PZN
CH 10	1,5 ml	40 cm	01/03/10/10/01/GHC	13851171
CH 12	2,0 ml	40 cm	01/03/10/12/02/GHC	13851188
CH 14	3,0 ml	40 cm	01/03/10/14/03/GHC	13851202
CH 16	5,0 ml	40 cm	01/03/10/16/05/GHC	13851219
CH 18	5,0 ml	40 cm	01/03/10/18/05/GHC	13851225

Packaging Unit: 10 pieces

Packaging: sterile

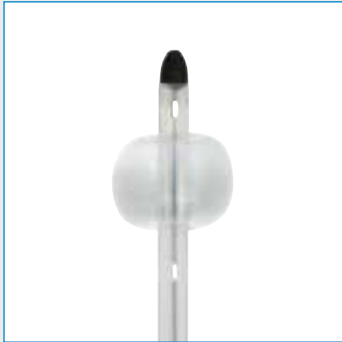

Balloon volume: 1,5 - 5,0 ml

CARE SET
Catheter insertion kit Basic and Clinic

Product specifications:

- ✓ The Basic Set comprises an inner and an outer package.
- ✓ The outer package of the basic set contains:
one sterile 10 ml syringe for unblocking the placed catheter.
- ✓ The items in the inner, sterilised package are used to place a new transurethral catheter.



Article	Article number	PZN
Basic	01/01/99/BS/GHC	01972188
Clinic	01/01/99/BSK/GHC	08845872

CARE SET
Catheter insertion kit Comfort and Comfort Plus

Product specifications:

- ✓ The Catheter insertion Kits comprises an inner and an outer package.
- ✓ The outer package of the basic set contains: one sterile 10 ml syringe for unblocking the inserted catheter a mucous membrane disinfectant Octenisept 15 ml.
- ✓ Additional included in Comfort Plus: lubricant Endosgel 6 ml.
- ✓ The items in the inner, sterilised package are used to place a new transurethral catheter.



Article	Article number	PZN
Comfort	01/01/99/KL/CO/GHC	14209815
Comfort Plus	01/01/99/KL/CP/GHC	14209821

Components of the catheter insertion kit are:

1 Tray



6 Swabs



1 Pair of tweezers



2 Underpads



1 Blocker syringe with 10% glycerine



3 Nitrile gloves



1 Fenestrated drape



Packaging Unit: 20 pieces
Packaging: sterile

Additional included in Basic set:



10 ml syringe

Pictures may differ from the original

Components of the catheter insertion kit are:

1 Tray



6 Swabs



1 Pair of tweezers



10 ml Syringe



1 Blocker syringe with 10% glycerine



3 Nitrile gloves



1 Fenestrated drape



2 Underpads



Mucous membrane disinfectant



Octenisept 15 ml

Additional included in Comfort Plus:



Lubricant Endosgel 6 ml

Packaging Unit: 35 pieces
Packaging: sterile

Pictures may differ from the original

CARE SET
Catheter insertion kit Comfort Eco

Product specifications:

- ✓ The Comfort Eco Set comprises an inner and an outer package.
- ✓ The outer package of the basic set contains: one sterile 10 ml syringe for unblocking the inserted catheter, mucous membrane disinfectant Octenisept 15 ml, lubricant LUBRICARE 6 ml.
- ✓ The items in the inner, sterilised package are used to place a new transurethral catheter.



Article	Article number	PZN
Comfort Eco	01/01/99/KL/CE/GHC	16134211

Additional equipment

Accessories for the catheter replacement

Next to our assortment covering catheters as well as drainage- and holding- systems, we also provide necessary accessories, to maintain a satisfactory supply of revulsive incontinence aids and ensuring an optimal replacement of a catheter, which all helps to maintain a high level of supply quality.

We want to highlight:

- The **blocking syringe containing a 10% glycerine solution, 10 ml** (Art.-Nr. 01/01/99/KB/10/GHC), ensures a steady volume in the catheter balloon during the wearing period, by preventing the liquid medium from diffusing out of the balloon. Thereby the balloon catheter stays fixated, and the risk of a dislocated catheter can almost be excluded.
- The **lubricant LUBRICARE, 6 ml** (Art.-Nr. 01/01/99/GG/001), adhere excellently onto the mucous membrane and has excellent lubrication properties, which ensures a smooth replacement of a catheter.

The below illustrated products are only an extract of our comprehensive assortment. We gladly provide more information upon request.




Blocker syringe with 10% glycerine, 10 ml



Lubricant LUBRICARE, 06 ml



Catheter plug

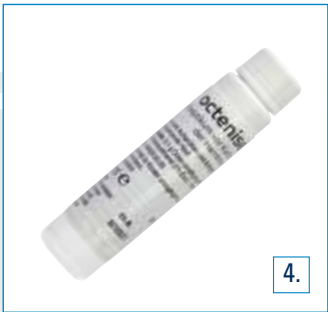

Octenisept, 15 ml

Components of the catheter insertion kit are:

1 Tray	6 Swabs	1 Pair of tweezers	10 ml Syringe	Mucous membrane
1 Blocker syringe with 10% glycerine	3 Nitrile gloves	1 Fenestrated drape	2 Underpads	Octenisept, 15 ml
				Lubricant
				LUBRICARE, 06 ml

Packaging Unit: 35 pieces
Packaging: sterile

Article	Article number	PZN	VE
1. Blocker syringe	01/01/99/KB/10/GHC	06718980	1VE = 200 pieces
2. LUBRICARE Lubricant	01/01/99/GG/06/001	13360220	1VE = 100 pieces
LUBRICARE Lubricant	01/01/99/GG/06/25	16938435	1VE = 4 x 25 pieces
Endosgel Lubricant	01/01/99/HG/EG/GHC	02437623	1VE = 250 pieces
3. Catheter plug	01/01/99/ZB/ST/GHC	14354919	1VE = 20 pieces
4. Octenisept	01/01/99/OCT/15/GHC	04804008	1VE = 100 pieces

Packaging: sterile

Pictures may differ from the original

CARE CATH_{Liv}

IC disposable catheter, hydrophilic, 17 cm

Ready for immediate use, Nelaton

Product specifications:

- ✓ Very good gliding properties thanks to hydrophilic coating.
- ✓ Ready for immediate use thanks to ready-to-use, hydrophilic coating.
- ✓ Atraumatic application due to rounded drainage eyes.
- ✓ Optimal length of 17 cm ensures complete bladder emptying.
- ✓ Non-contact application and good guideability due to the correct rigidity of the catheter.
- ✓ Good grip of the handle allows good handling, even with limited dexterity.
- ✓ Easy connection of urine collection bag and/or instillation solution via the handle cone.



Atraumatic, rounded drainage eyes



Grip of the handle

CARE CATH_{Liv}

Size	Length	Item	number
CH 10	17 cm	LIV-H1017	18890885
CH 12	17 cm	LIV-H1217	18890916
CH 14	17 cm	LIV-H1417	18890922

Packaging unit: 30 pieces
Packaging: sterile

Illustrations may differ from the original

CARE CATH_{Liv}

IC disposable catheter, hydrophilic, 41 cm

Ready for immediate use, Nelaton

Product specifications:

- ✓ Very good gliding properties thanks to hydrophilic coating.
- ✓ Ready for immediate use thanks to ready-to-use, hydrophilic coating.
- ✓ Atraumatic application due to rounded drainage eyes.
- ✓ Non-contact application and good guideability due to the correct rigidity of the catheter.
- ✓ Good grip of the handle allows good handling, even with limited dexterity.
- ✓ Easy connection of urine collection bag and/or instillation solution via the handle cone.



Atraumatic, rounded drainage eyes



Grip of the handle

CARE CATH_{Liv}

Size	Length	Item	number
CH 12	41 cm	LEIV-H1241	18890833
CH 14	41 cm	LEIV-H1441	18890856
CH 16	41 cm	LEIV-H1641	18890862

Packaging unit: 30 pieces
Packaging: sterile

Pictures may differ from the original

CARE CON comfort | Typ I
Male External Catheter (MEC) 100% silicone
Self-adhesive

Product specifications:

- ✓ The CARE CON comfort captivates through its comfortable silicone material and good hold.
- ✓ The material thickness of the urinary sheath has been developed to cling softly and flexibly whilst ensuring a secure hold, which provides high safety.
- ✓ The stable connecting tube prevents the MEC from bending and ensures a secure connection to the urine bag.
- ✓ Small steps between the sizes ensure perfect fit.
- ✓ The medical grade silicone maintains the breathing activity of the skin.
- ✓ Due to the use of transparent materials, skin alterations or irritations can be observed at any time.
- ✓ With a shaft length of 8,0 cm and an adhesive area of 4,5 cm, the CARE CON comfort meets the requirements of most MEC-users.



CARE CON special | Type II
Male External Catheter (MEC) 100% silicone
Self-adhesive

Product specifications:

- ✓ The CARE CON special captivates through its comfortable silicone material and good hold, despite the short shaft length.
- ✓ The material thickness of the urinary sheath has been developed to cling softly and flexibly whilst ensuring a secure hold, which provides high safety.
- ✓ The stable connecting tube prevents the MEC from bending and ensures a secure connection to the urine bag.
- ✓ Small steps between the sizes ensure perfect fit.
- ✓ With a shaft length of 5,0 cm and an adhesive area of 4,0 cm, the CARE CON special meets the special requirements of users with a retracted penis.
- ✓ The medical grade silicone maintains the breathing activity of the skin.
- ✓ Due to the use of transparent materials, skin alterations or irritations can be observed at any time.



Diameter	Shaft / Adhesive surface	Article number	PZN
25 mm	8,0 cm / 4,5 cm	01/01/28/01/25/GHC	17520554
28 mm	8,0 cm / 4,5 cm	01/01/28/01/28/GHC	17520583
30 mm	8,0 cm / 4,5 cm	01/01/28/01/30/GHC	17520620
32 mm	8,0 cm / 4,5 cm	01/01/28/01/32/GHC	17520637
36 mm	8,0 cm / 4,5 cm	01/01/28/01/36/GHC	17520643
41 mm	8,0 cm / 4,5 cm	01/01/28/01/41/GHC	17520666

Packaging Unit: 30 pieces
Packaging: non-sterile

Pictures may differ from the original

Diameter	Shaft / Adhesive surface	Article number	PZN
25 mm	5,0 cm / 4,0 cm	01/01/28/02/25/GHC	17520672
28 mm	5,0 cm / 4,0 cm	01/01/28/02/28/GHC	17520689
30 mm	5,0 cm / 4,0 cm	01/01/28/02/30/GHC	17520703
32 mm	5,0 cm / 4,0 cm	01/01/28/02/32/GHC	17520726
36 mm	5,0 cm / 4,0 cm	01/01/28/02/36/GHC	17520732
41 mm	5,0 cm / 4,0 cm	01/01/28/02/41/GHC	17520749

Packaging Unit: 30 pieces
Packaging: non-sterile

Pictures may differ from the original

CARE CON active | Typ III

Male External Catheter (MEC) 100% silicone

Self-adhesive

Product specifications:

- ✓ The CARE CON active captivates through its comfortable silicone material and good hold, with an extra-long adhesive area.
- ✓ The material thickness of the urinary sheath has been developed to cling softly and flexibly whilst ensuring a secure hold, which provides high safety.
- ✓ The stable connecting tube prevents the MEC from bending and ensures a secure connection to the urine bag.
- ✓ Small steps between the sizes ensure perfect fit.
- ✓ With a shaft length of 8,0 cm and an adhesive area of 7,0 cm, the CARE CON active meets the requirements of sporty and active MEC-users.
- ✓ The medical grade silicone maintains the breathing activity of the skin.
- ✓ Due to the use of transparent materials, skin alterations or irritations can be observed at any time.



Diameter	Shaft / Adhesive surface	Article number	PZN
25 mm	8,0 cm / 7,0 cm	01/01/28/03/25/GHC	17520755
28 mm	8,0 cm / 7,0 cm	01/01/28/03/28/GHC	17520761
30 mm	8,0 cm / 7,0 cm	01/01/28/03/30/GHC	17520784
32 mm	8,0 cm / 7,0 cm	01/01/28/03/32/GHC	17520790
36 mm	8,0 cm / 7,0 cm	01/01/28/03/36/GHC	17520809
41 mm	8,0 cm / 7,0 cm	01/01/28/03/41/GHC	17520815

Packaging Unit: 30 pieces

Packaging: non-sterile

Pictures may differ from the original

CARE FLOW UB 4000TUR closed TUR irrigation bag system

with pump ball for controlled post-operative bladder irrigation

Capacity: 4000 ml

Product specifications:

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder with fixing cord.
- ✓ Bacteria- and fluid-tight venting filters ensure permanent pressure compensation and provide for optimal drainage in the system.
- ✓ Due to controlled irrigation of the bladder with the pump ball integrated in the feeding tube, post-operative tamponade of the bladder can be avoided.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.



Capacity	Tube length	Article number	PZN
4000 ml	120 cm	04/02/4000/TUR/GHC	06189004



Pump ball

The specific interaction of the pump ball and the two tube clamps allows to flush out blood clots in the bladder and tube system.



Drip chamber

An anti-reflux valve divides the drip chamber and reliably prevents a backflow of urine and thereby obviates the risk against an ascending infection.



Profiled tube

The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



Drain tap in the doghouse

The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

Packaging Unit: 10 pieces

Packaging: sterile

Pictures may differ from the original

CARE FLOW UB 4000

Closed urine drainage system

Capacity: 4000 ml

Product specifications:

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder with fixing cord.
- ✓ Bacteria- and fluid-tight venting filters ensure permanent pressure compensation and provide for optimal drainage in the system.
- ✓ Clear graduation markings allow for reliable volume control.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.
- ✓ T-drain tap.
- ✓ Kink resistant tube.



Capacity	Tube length	Article number	PZN
4000 ml	120 cm	04/02/4000/120/GHC	06718939

CARE FLOW UB 2000

Closed urine drainage system

Capacity: 2000 ml

Product specifications:

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder with fixing cord.
- ✓ Bacteria- and fluid-tight venting filters ensure permanent pressure compensation and provide for optimal drainage in the system.
- ✓ Clear graduation markings allow for reliable volume control.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.
- ✓ T-drain tap.
- ✓ Kink resistant tube.



Capacity	Tube length	Article number	PZN
2000 ml	120 cm	04/02/2000/120/GHC	06718922
2000 ml	90 cm	04/02/2000/90/GHC	09151597

Packaging Unit: 20 pieces
Packaging: sterile

The **CARE FLOW BAG** is a cover for the **CARE FLOW UB 2000** serving as a discreet way to store the urinary bedside bag.

Product specifications:

- ✓ Thanks to the convenient Velcro fastening the cover can be affixed to the urinary bag quickly and easily. A second Velcro fastening on the front allows for a clean and uncomplicated extraction and reinsertion of the drain tap.
- ✓ The cover is made from 100% polypropylene (PP) and is thus resistant against organic acids.
- ✓ The cover is light grey, washable at 30° and hence reusable.

Article number	PZN
04/99/01/945/GHC	09764331

Packaging Unit: 10 pieces
Packaging: non-sterile



Staggered connector
 The staggered connector provides a reliable connection with the catheter. A needle-free sampling of urine is possible at any time, due to the needle-free sample port, which prevents needlestick injuries.



Profiled tube
 The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



Drip chamber
 An anti-reflux valve divides the drip chamber and reliably prevents a backflow of urine and thereby obviates the risk against an ascending infection.



Drain tap in the doghouse
 The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

Packaging Unit: 20 pieces
Packaging: sterile

Pictures may differ from the original

CARE FLOW UB 2000A Disposable bedside bag

Sterile bedside bag with integrated twin holder

Capacity: 2000 ml

Product specifications:

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.
- ✓ Clear graduation markings allow for reliable volume control.
- ✓ A slit ensures the safe and hygienic positioning of the dRAINTap.
- ✓ T-drain tap.
- ✓ Kink resistant tube.



Capacity	Tube length	Article number	PZN
2000 ml	110 cm	04/04/2000/110/GHC	08871289

CARE FLOW UB 2000E, UB 2000US

Sterile / non-sterile single-use bedside bag with holder

Capacity: 2000 ml

Product specifications:

- ✓ The system is equipped with a universal staggered cone and can be connected to all common catheters.
- ✓ The product must be changed daily.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.
- ✓ Use the twin holder shown below to attach the bag to many different types of bed frames.
- ✓ The urine bag is emptied safely by means of the T-valve.



Article	Capacity	Sterility	Tube length	Article number	PZN
UB 2000E	2000 ml	sterile	90 cm	04/04/2000/90/GHC	16047274
UB 2000US	2000 ml	non-sterile	90 cm	04/06/2000/90/GHC	12418549
UB 2000US	2000 ml	non-sterile	130 cm	04/06/2000/130/GHC	08845487

UB 2000E: **Packaging Unit:** 100 pieces, **Packaging:** sterile
 UB 2000US: **Packaging Unit:** 100 pieces, **Packaging:** non-sterile **Sub-packaging:** each 10 pieces

Holder for UB 2000E and UB 2000US

- ✓ Plastic hanger for bedside bags with twin holder. Equipped with two pins taking up single-use bedside bags hanger holes.

Article number	PZN
04/99/BHL/PVC/GHC	06718951

Packaging Unit: 100 pieces **Sub-packaging:** 10 pieces each
Packaging: non-sterile



46



Staggered connector
 The staggered connector provides a reliable connection with the catheter. A needle-free sampling of urine is possible at any time, due to the needle-free sample port, which prevents needlestick injuries.



Profiled tube
 The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



Venting filter
 A bacteria and liquid-tight venting filters ensure permanent pressure equalisation. This ensures an optimal drainage of the system.



Drain tap
 The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

Packaging Unit: 50 pieces
Packaging: sterile

Pictures may differ from the original

47

Pictures may differ from the original

CARE FLOW UB 2000USo, UB 1500USo, SB 1500USo

Non-sterile single-use bedside bag/secretion bag without drain tap

Capacity: 2000 ml, 1500 ml

Product specifications:

- ✓ The system is equipped with a universal staggered cone and can be connected to all common catheters.
- ✓ The product must be changed daily.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria. (except the SB 1500USo)
- ✓ Use the twin holder shown below to attach the bag to many different types of bed frames.



Article	Capacity	Tube length	Anti-reflux valve	Article number	PZN
UB 2000USo	2000 ml	90 cm	yes	04/07/2000/90/GHC	08845501
UB 1500USo	1500 ml	90 cm	yes	04/07/1500/90/01	16839791
SB 1500USo	1500 ml	90 cm	no	04/07/1500/90/00	16839785

CARE FLOW UB 1000 Closed urine drainage systems

Leg bag systems for day & night

Capacity: 1000 ml

Product specifications:

- ✓ Closed urine drainage system
- ✓ Bacteria- and fluid-tight venting filters provide for permanent pressure compensation and optimal drainage in the system.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.
- ✓ T-drain tap.
- ✓ Needle free sampling port.
- ✓ Clear gradvation marks.



Capacity	Tube length	Article number	PZN
1000 ml	30 cm	04/03/1000/30/GHC	06718885
1000 ml	60 cm	04/03/1000/60/GHC	06718891
1000 ml	90 cm	04/03/1000/90/GHC	06718916

48

UB 2000USo:

Packaging Unit: 150 pieces, Packaging: non-sterile

Sub-packaging: each 10 pieces

UB 1500USo:

Packaging Unit: 200 pieces, Packaging: non-sterile

Sub-packaging: each 10 pieces

SB 1500USo:

Packaging Unit: 200 pieces, Packaging: non-sterile

Sub-packaging: each 10 pieces

Holder for UB 2000USo, UB 1500USo and SB 1500USo

- ✓ Plastic hanger for bedside bags with twin holder. Equipped with two pins taking up single-use bedside bags hanger holes.



Article number	PZN
04/99/BHL/PVC/GHC	06718951

Packaging Unit: 100 pieces Sub-packaging: pieces each

Packaging: non-sterile

Pictures may differ from the original

49



Stepped connector

The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



Fleece coating

The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



Drip chamber

An anti-reflux valve in the drip chamber provides for reliable prevention of urine reflux and highly efficient protection against ascending migration of bacteria.



Drain tap

The special design of the piston in the drain tap ensures safe and easy operation and prevents leak-ages caused by improper use of the piston.

Packaging Unit: 30 pieces

Packaging: sterile

Pictures may differ from the original

Fixation for leg bags on page 51-56.

CARE FLOW UB 750 + 750US sterile / non-sterile

1-way chamber leg bag systems for mobile users

Capacity: 750 ml

Product specifications:

- ✓ The leg bag with 10 cm and 30 cm tube can conveniently be worn on the thigh.
- ✓ The leg bag with 50 cm tube comes with an additional stepped connector in the package. So the tube can be cut to the individual length and the leg bag can conveniently be worn on the thigh or lower leg.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.

Capacity	Tube length	Article number	PZN	Sterility
750 ml	10 cm	04/01/750/10/GHC	06718862	sterile
750 ml	30 cm	04/01/750/30/GHC	09772046	sterile
750 ml	50 cm	04/01/750/50/GHC	06718879	sterile

Packaging Unit: 10 cm: 75 pieces 30 cm: 60 pieces, 50 cm: 50 pieces

Capacity	Tube length	Article number	PZN	Sterility
750 ml	10 cm	04/01/750US/10/GHC	01972716	non-sterile
750 ml	30 cm	04/01/750US/30/GHC	12341427	non-sterile
750 ml	50 cm	04/01/750US/50/GHC	00226543	non-sterile

Packaging Unit: 10 cm: 120 pieces 30 cm: 100 pieces, 50 cm: 100 pieces, Sub-Packaging: each 10 pieces



CARE FLOW UB 750F + UB 750UF sterile / non-sterile

1-way chamber leg bag systems with flip valve for mobile users

Capacity: 750 ml

Product specifications:

- ✓ Especially users with limited hand dexterity or reduced power can safely and hygienically empty this leg bag, due to our newly developed flip valve.
- ✓ The leg bag with 50 cm tube comes with an additional stepped connector in the package. So the tube can be cut to the individual length and the leg bag can conveniently be worn on the thigh or lower leg.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.

Capacity	Tube length	Article number	PZN	Sterility
750 ml	50 cm	04/01/750F/50/GHC	12371500	sterile

Packaging Unit: 50 cm: 50 pieces

Capacity	Tube length	Article number	PZN	Sterility
750 ml	50 cm	04/01/750UF/50/GHC	12371517	non-sterile

Packaging Unit: 50 cm: 100 pieces

Packaging: each 10 pieces



50



Stepped connector

The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



Doghouse

A doghouse on the bag ensures safe and hygienic positioning of the drain tap.



Fleece coating

The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



Drain tap

The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

Packaging: sterile / non-sterile

Fixation for leg bags on page 51-56.

Pictures may differ from the original

51



Stepped connector

The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



Doghouse

A doghouse on the bag ensures safe and hygienic positioning of the drain tap.



Fleece coating

The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



Drainage tap

The ergonomically formed flat flip valve of the drainage tube can be used easily and intuitively. To drain the leg bag the lever of the valve is flipped downwards by 180°.

Packaging: sterile / non-sterile

Fixation for leg bags on page 51-56.

Pictures may differ from the original

CARE FLOW UB 500 + UB 500US sterile / non-sterile

3-way chamber leg bag system for mobile users

Capacity: 500 ml

Product specifications:

- ✓ The leg bag with 10 cm and 30 cm tube can conveniently be worn on the thigh.
- ✓ The leg bag with 50 cm tube comes with an additional stepped connector in the package. So the tube can be cut to the individual length and the leg bag can conveniently be worn on the thigh or lower leg.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.

Capacity	Tube length	Article number	PZN	Sterility
500 ml	10 cm	04/01/500/10/GHC	09084990	steril
500 ml	30 cm	04/01/500/30/GHC	12341203	steril
500 ml	50 cm	04/01/500/50/GHC	09085009	steril

Packaging unit: 10 cm: 75 pieces, 30 cm: 60 pieces, 50 cm: 50 pieces

Capacity	Tube length	Article number	PZN	Sterility
500 ml	10 cm	04/01/500US/10/GHC	09085015	non-sterile
500 ml	30 cm	04/01/500US/30/GHC	12341410	non-sterile
500 ml	50 cm	04/01/500US/50/GHC	09085021	non-sterile

Packaging unit: 10 cm: 120 pieces, 30 cm: 100 pieces, 50 cm: 100 pieces, Sub-Packaging: each 10 pieces



CARE FIX

Universal Velcro strap

Length can be individually adjusted

Product specifications:

- ✓ Universal leg bag holding band with two attached Velcro straps to go through the upper eyelets of the leg bag and a special sewn on fixation for the tube of the leg bag system.
- ✓ The elastic material ensures that disturbances of the blood circulation in the thigh and lower leg are prevented as far as possible.
- ✓ The Velcro strap is 70 cm long and can be cut to the individual length. So it can be used on the thigh and lower leg.
- ✓ This universal Velcro strap can be used for the following urine leg bag systems: UB1000, UB750/UB750US, UB750F/UB750UF, UB500/UB500US.



Stepped connector

The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



3-chambers

Due to the 3-chambers will the urin spread more regular through the system. In that way is it possible to wear the filled bag even more discrete beneath the clothing.



Fleece coating

The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



Drainage tap

The ergonomically formed flat flip valve of the drainage tube can be used easily and intuitively. To drain the leg bag the lever of the valve is flipped downwards by 180°.

Packaging: sterile / non-sterile

Fixation for leg bags on page 51-56.

Pictures may differ from the original

- 60°
- Do not put in the dryer
- Do not iron

Lenght	Article number	PZN
70 cm	05/01/BKL/Uni/GHC	06718968

Packaging Unit: 30 pieces
Packaging: sterile



Pictures may differ from the original

CARE FIX

Velcro strap premium set

Length can be individually adjusted

Product specifications:

- ✓ Sets consists of:
an extra wide holding band with two sewn-on
Velcro straps to fixate the upper part of the leg bag.
- ✓ Both sets are equipped sewn-on loop to fixate the
leg bag tube.
- ✓ The second Velcro Strap ensures a safe fixation of
the lower part of the leg bag.
- ✓ Silicone strips on the inner side of the holding band
ensures a good grip on the leg.
- ✓ The elastic material prevents a disruption of the blood
circulation on the thigh or lower leg.
- ✓ A 70 cm long and individually adjustable Velcro strap can
be worn on both the thigh as well as the lower leg.
- ✓ The attached loop for an easier opening of the holding
band.

Lenght	Article number	PZN
70 cm	05/01/BKL/PRE/SET	09241099

Packaging Unit: 10 pieces
Packaging: non-sterile



CARE FIX

Leg pants with integrated leg bag pocket

Product specifications:

- ✓ Soft and ribbed lag pants with integrated
leg bag pocket.
- ✓ Comfortable and safe fixation on the High
as well as the lower leg.
- ✓ Due to the special design of the drain tap on
the leg bag, emptying the bag is very simple.
- ✓ The extra-wide lateral Velcro fastener makes
it very easy for the user to put on the leg pan
(Art.-Nr. 05/02/10).
- ✓ The CARE FIX leg pants can be used for the
following urine leg bag systems:
UB1000, UB750/UB750US, UB750F/UB750UF,
UB500/UB500US.

Size	Article number	PZN
XXL	05/02/10/OS/XXL/GHC	08880957
XL	05/02/10/OS/XL/GHC	08871272
L	05/02/10/OS/L/GHC	08876513
M	05/02/10/OS/M/GHC	08845748
S	05/02/10/US/S/GHC	07126649

Measurement: XXL = 60 - 85 cm, XL = 50 - 60 cm,
L = 40 - 60 cm, M = 35 - 45 cm, S = 25 - 40 cm

Packaging Unit: 10 pieces
Packaging: non-sterile

- 60°
- Do not put in the dryer
- Do not iron

Measure in the middle of the upper or lower leg.

Size	Article number	PZN
XL	05/02/11/OS/XL/GHC	12418555
L	05/02/11/OS/L/GHC	07126017
M	05/02/11/OS/M/GHC	07126023
S	05/02/11/US/S/GHC	07126313

Measurement: XL = 65 - 95 cm, L = 55 - 80 cm,
M = 45 - 65 cm, S = 30 - 45 cm

Packaging Unit: 10 pieces
Packaging: non-sterile



54

- ✓ This Velcro strap premium set can be used for the
following urine leg bag systems: UB1000, UB750/
UB750US, UB750F/UB750UF, UB500/UB500US.
- ✓ The premium set with buttons (Art.-Nr. 05/01/BKL/PRE/03)
comes with additional four buttons to fixate the upper
part of the leg bag.

- 60°
- Do not put in the dryer
- Do not iron

Lenght	Article number	PZN
70 cm	05/01/BKL/PRE/03	12357931

Packaging Unit: 10 pieces
Packaging: non-sterile



55

Pictures may differ from the original

Pictures may differ from the original

CARE FIX
Thigh cuffs sizes: XL, L, M

Product specifications:

- ✓ The fabric leg bag holder made from a washable cotton blend has an integrated leg bag pocket.
- ✓ A silicon strip on the inner side ensures a good grip.
- ✓ The fabric leg bag holder can be worn on both the right as well as the left thigh.
- ✓ The bags are positioned into the integrated pocket with a Velcro fastener; including an opening for the drain tap.
- ✓ The CARE FIX cuffs can be used for the following urine leg bag systems: UB1000, UB750/UB750US, UB750F/UB750UF, UB500/UB500US.



CARE FIX
Lower leg cuffs sizes: L, M, S

Product specifications:

- ✓ The fabric lower leg bag holder made from a washable cotton blend has an integrated leg bag pocket.
- ✓ A silicon strip on the inner side ensures a good grip.
- ✓ The fabric lower leg bag holder can be worn on both the right as well as the left lower leg.
- ✓ Leg bags are positioned into the integrated pocket with a Velcro fastener; including an opening for the drain tap.
- ✓ The CARE FIX cuffs can be used for the following urine leg bag systems: UB1000, UB750/UB750US, UB750F/UB750UF, UB500/UB500US.



56

- 60°
- Do not put in the dryer
- Do not iron

It is measured in the middle of the thigh.

Leg bag holder for the thigh

Size	Article number	PZN
XL	05/02/12/XL/GHC	09080408
L	05/02/12/L/GHC	09080383
M	05/02/12/M/GHC	09080377

Measurement: XL = 61 cm, L = 58 cm, M = 54 cm

Packaging Unit: 10 pieces

Packaging: non-sterile



Pictures may differ from the original

57

- 60°
- Do not put in the dryer
- Do not iron

It is measured in the middle of the lower leg.

Lower Leg Cuffs

Size	Article number	PZN
L	05/02/13/L/GHC	09080437
M	05/02/13/M/GHC	09080420
S	05/02/13/S/GHC	09080414

Measurement: L = 37 cm, M = 35 cm, S = 32 cm

Packaging Unit: 10 pieces

Packaging: non-sterile



Pictures may differ from the original

CARE FIX
Premium one leg pants sizes: XL, L, M

Product specifications:

- ✓ The premium one leg pants are made from a washable cotton blend has an integrated leg bag pocket.
- ✓ Individually adjustable belt with silicone strips for a safe fixation.
- ✓ The premium one leg pant can be worn on both the right as well as the left thigh.
- ✓ Leg bags are positioned into the integrated pocket with a Velcro fastener; including an opening for the drain tap.
- ✓ Easy emptying of the leg bags due to the comfortable positioning of the drain tap.
- ✓ The leg pants can be used for the following urine leg bag systems: UB1000, UB750/UB750US, UB750F/UB750UF, UB500/UB500US.



- 60°
- Do not put in the dryer
- Do not iron

Measure in the middle of the upper or lower leg.

Size	Article number	PZN
XL	05/03/10/XL/GHC	09080472
L	05/03/10/L/GHC	09080466
M	05/03/10/M/GHC	09080443

Measurement: XL = 98 - 120 cm, L = 90 - 100 cm, M = 85 - 95 cm

Packaging Unit: 10 pieces

Packaging: non-sterile



Pictures may differ from the original

CARE FLOW und CARE VENT catheter valve
For an independent and controlled emptying of the urinary bladder

Product specifications:

- ✓ The optimal solution for mobile, independent users wearing a transurethral or suprapubic catheter.
- ✓ The catheter valve continuously trains the bladder muscles so that the bladder function remains completely intact. The emptying of the urinary bladder follows the user's personal rhythm. This sustainably enhances the quality of life.
- ✓ During the night the valve can be connected to a urine bag system by means of a firmly attached silicone tube.
- ✓ By flipping the lever of the flip valve by 180° the CARE FLOW catheter valve is completely opened as well as fixed in a safe position. This prevents an unintentional closing of the valve during the night.



Arrows indicate open / closed position



Closed catheter valve



Open catheter valve, flipped by 180°

CARE FLOW catheter valve

Weight	Article number	PZN
3g	03/02/10/00/GHC	12378382

Packaging Unit: 5 pieces

Packaging: sterile

- ✓ The special construction prevents the after-dripping of urine, due to the which ensures a hygienical operation and utmost discretion.
- ✓ A colour coding of the pivot and slide mechanism makes it very easy to use single handed.

CARE VENT catheter valve

Weight	Article number	PZN
6g	03/02/10/01/GHC	01914676

Packaging Unit: 10 pieces

Packaging: non-sterile



Pictures may differ from the original

CARE DRESS
Suprapubic dressing set

Product specifications:

- ✓ Suprapubic dressing set for optimal care of the suprapubic exit site.
- ✓ The included smaller plaster serves to temporarily fix the suprapubic catheter to the abdominal wall, thus preventing the catheter from shifting.



CARE DRESS
Abdominal PEG dressing set

Product specifications:

- ✓ Abdominal PEG dressing set for optimal care of the PEG exit site.



Included in this set:

- 1 x elastic plaster 10 x 10 cm
- 1 x elastic plaster 5 x 3 cm
- 3 x non-woven swabs 5 x 5 cm, 6-ply
- 2 x non-woven split swabs 5 x 5 cm, 6-ply
- Packed in a peel-apart pouch -

Article number	PZN
08/02/VS/GHC	06718974

Packaging Unit: 25 sets
Packaging: sterile

Pictures may differ from the original

Included in this set:

- 1 x elastic plaster 12 x 15 cm
- 3 x non-woven swabs 7,5 x 7,5 cm, 8 ply
- 2 x non-woven split swabs 7,5 x 7,5 cm, 4 ply
- Packed in a peel-apart pouch -

Article number	PZN
08/02/PEG/GHC	05461059

Packaging Unit: 15 sets
Packaging: sterile

Pictures may differ from the original



Visit our website www.ghc-medical.de.

There you will find detailed information about our company as well as the entire product portfolio.

Our sales team is looking forward to your call:

Mon.-Thurs.: 08:00 – 17:00
Fri.: 08:00 – 14:00

Contact:

GHC German Health Care GmbH
Birkstück 15
24999 Wees/Flensburg

Tel.: +49 (0) 46 31 / 57 200-0
Fax: +49 (0) 46 31 / 57 200-19

info@ghc-medical.de
www.ghc-medical.de