


# BLUECATH

## Suprapubic *integral* 100% silicone catheter

### 2-way integral balloon catheter with plug

#### Product specifications:

- ✓ The catheter is equipped with an integral balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ Safe fixation of the catheter, due to the integral balloon.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short, central open tip.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.

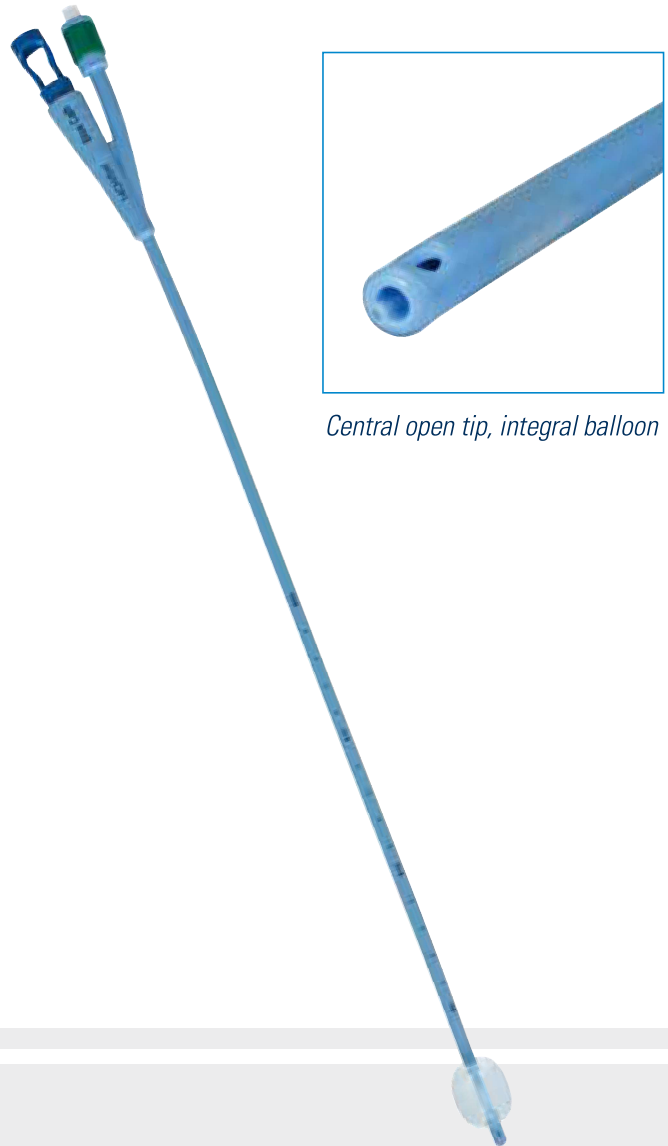

Central open tip, integral balloon

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/02/22/12/05/GHC	11851037
CH 14	05 ml	40 cm	01/02/22/14/05/GHC	11851043
CH 16	10 ml	40 cm	01/02/22/16/10/GHC	11851072
CH 18	10 ml	40 cm	01/02/22/18/10/GHC	11851089
CH 20	10 ml	40 cm	01/02/22/20/10/GHC	11851095
CH 22	10 ml	40 cm	01/02/22/22/10/GHC	11851126
CH 24	10 ml	40 cm	01/02/22/24/10/GHC	11851132

**Packaging Unit:** 10 pieces

**Packaging:** sterile



Balloon volume: 05-10 ml

Pictures may differ from the original