

CARE FLOW UB 1000 Closed urine drainage systems

Leg bag systems for day & night

Capacity: 1000 ml

Product specifications:

- ✓ For a safe and comfortable fixation of the leg bags during the day various high quality products from the product family CARE FIX are available.
- ✓ A hanger for bedside bags (which can be ordered separately) in combination with the fixing cord attached to the system ensures the safe and proper fixation to various bed frames during the night.
- ✓ Bacteria- and fluid-tight venting filters provide for permanent pressure compensation and optimal drainage in the system.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.




Capacity	Tube length	Article number	PZN
1000 ml	30 cm	04/03/1000/30/GHC	06718885
1000 ml	60 cm	04/03/1000/60/GHC	06718891
1000 ml	90 cm	04/03/1000/90/GHC	06718916



Staggered connector

The staggered connector provides a reliable connection with the catheter. A needle-free sampling of urine is possible at any time, due to the needle-free sample port, which prevents needlestick injuries.

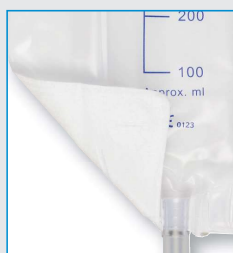

Fleece coating

The backside of the leg bag is equipped with a comfortable and sweat-absorbing fleece, in order to prevent skin irritations and it considerably enhances user comfort.



Drip chamber

An anti-reflux valve divides the drip chamber and reliably prevents a backflow of urine and thereby obviates the risk against an ascending infection.



Drain tap

The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

Packaging Unit: 30 pieces

Packaging: sterile

Fixation for leg bags on page 46-51.

Pictures may differ from the original