

# CARE FLOW UMS 2500

## Closed Urine measurement system, incl. CAREFLOW UB 2000UMS changing bag, 2000 ml

### Product specifications:

- ✓ The CAREFLOW UMS 2500 is perfect for precise measurement and collection of discharged urine.
- ✓ It consists of a measuring chamber with a volume of 500 ml as well as a removable urine bag with a capacity of 2000 ml.
- ✓ The twin holder with additional attachment cord as well as two-sided suspensions allow one secure attachment to a variety of bed throws.
- ✓ The bacteria and liquid-tight ventilation filters ensure permanent pressure equalization.
- ✓ The T-Tap drainage tap of the collection bag is one-handed and safe to use and allows a fast emptying.
- ✓ The needle-free sampling point prevents needlestick injuries.



### Urine measurement system

Capacity	Tube length	Article number	PZN
2500 ml	120 cm	04/08/2000/120/GHC	14323971
2500 ml	150 cm	04/08/2000/150/GHC	15423172

### Changing bag

Capacity	Article number	PZN
2000 ml	04/09/2000/GHC	14323965



#### Measuring chamber

The measuring chamber has a precise and clear graduation and is divided in a fine (1 ml) and coarse (10 ml) measuring unit.




#### Screw cap

With the screw cap of the measuring chamber a single-handed and a cross-contamination-free emptying of the collection bag is assured.



#### Attachment bag

The detachable bag is safely connected to the measuring chamber but can also easily be changed thanks to the bayonet lock and the two mountings where the bag is attached.

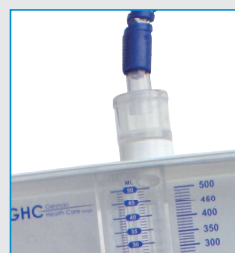

#### Drip chamber

An anti-reflux valve divides the drip chamber and reliably prevents a backflow of urine and thereby obviates the risk against an ascending infection.

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original